DSG Connect App Tutorial

English.....	2
Français.....	46
Deutsch.....	90
Português.....	134
Español.....	178

DSG Connect App Tutorial

English

1. What is DSG?.....	3-6	6. Parameters	33-41
2. The Interface Global Overview.....	7	DSG Device Connection.....	33
3. DSG Connect Device.....	8-9	Simple & Advanced Mode.....	34
4. Simple Mode.....	10-20	Volume.....	35
Signal.....	10-11	Label.....	36
Empty Fields - Notes.....	12-13	Signal.....	37
Reference lines.....	14-17	Languages.....	38
Screenshots.....	18-19	IFU.....	39
Pause.....	20	About Us.....	40
5. Advanced Mode	21-32	Exit App.....	41
Pedicle Selection.....	21-24	7. System Notifications.....	42-44
Recording.....	25-26	Battery.....	42
Label.....	27-29	Connection.....	43
Export-report (data review).....	30-32	Memory Capacity.....	44

What is DSG?

Dynamic Surgical Guidance (DSG[®]) Technology

The DSG Technology is based on the principle that the electrical conductivity of tissue varies by the type of tissue.

Electrical conductivity is the property of a material to conduct an electrical current. e.g.:

- Cortical bone has low conductivity.
- Cancellous bone has medium conductivity.
- Periosteum and blood have high conductivity.

The varying conductivity, is a parameter that can be very helpful in differentiating between various tissue types while drilling a pedicle prior to screw placement.

What is DSG?


How does the DSG Connect device work?

A bipolar sensor is embedded at the tip of each DSG Connect device probe. A current is pulsed through two electrodes in the probe. Based on the local conductivity, the sensor can accurately inform the surgeon of the type of tissue the tip is in. This in turn can alert the surgeon prior to an imminent cortical breach during pedicle preparation. The surgeon is informed in real time of the tissue type by changes in the pitch and cadence of an audio signal and a flashing LED light.

The sound signal will change depending whether the surgeon is in the cancellous bone, cortical bone or blood.


What is DSG?

Tip in cancellous bone

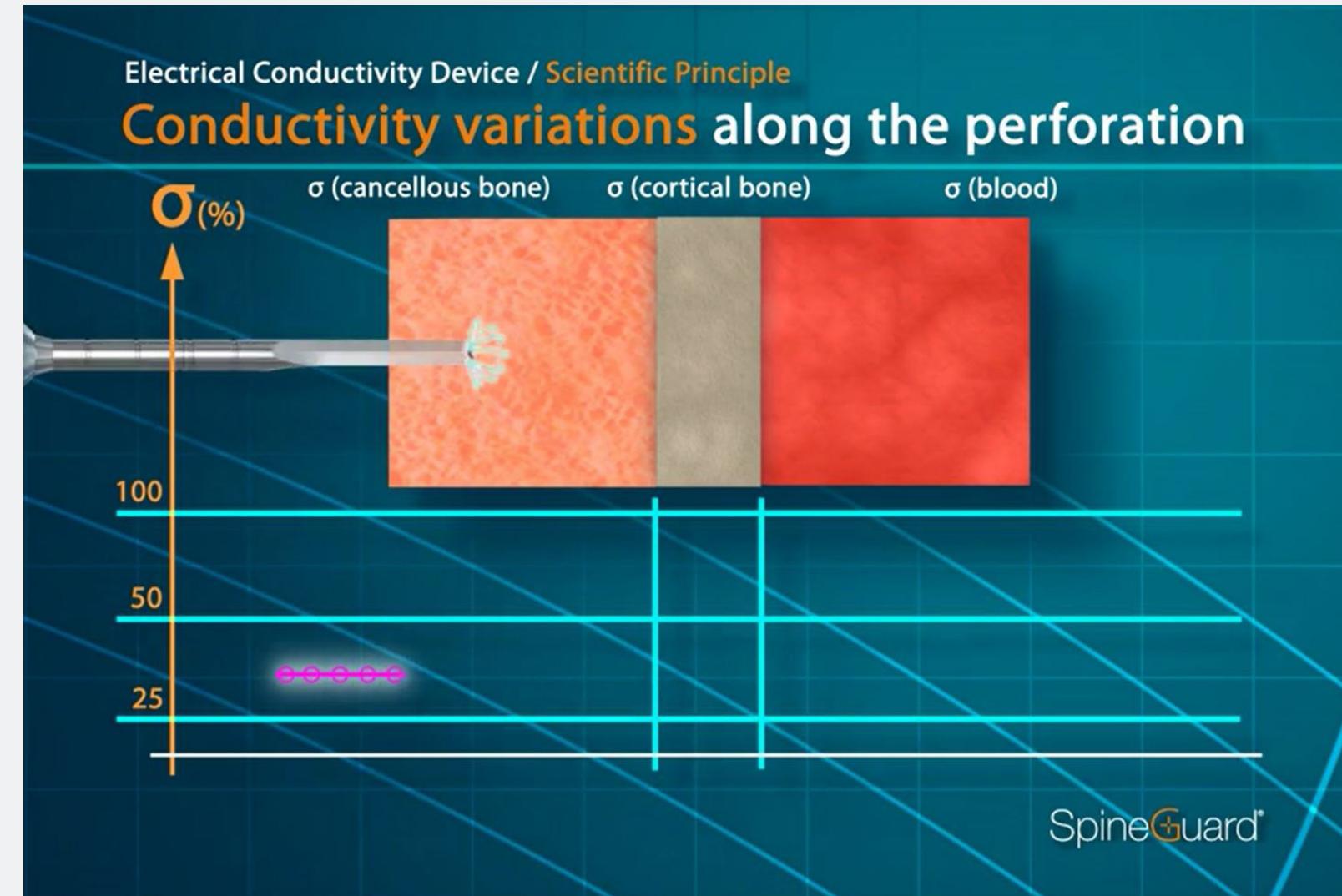
Medium pitch, medium cadence

Tip in cortical bone

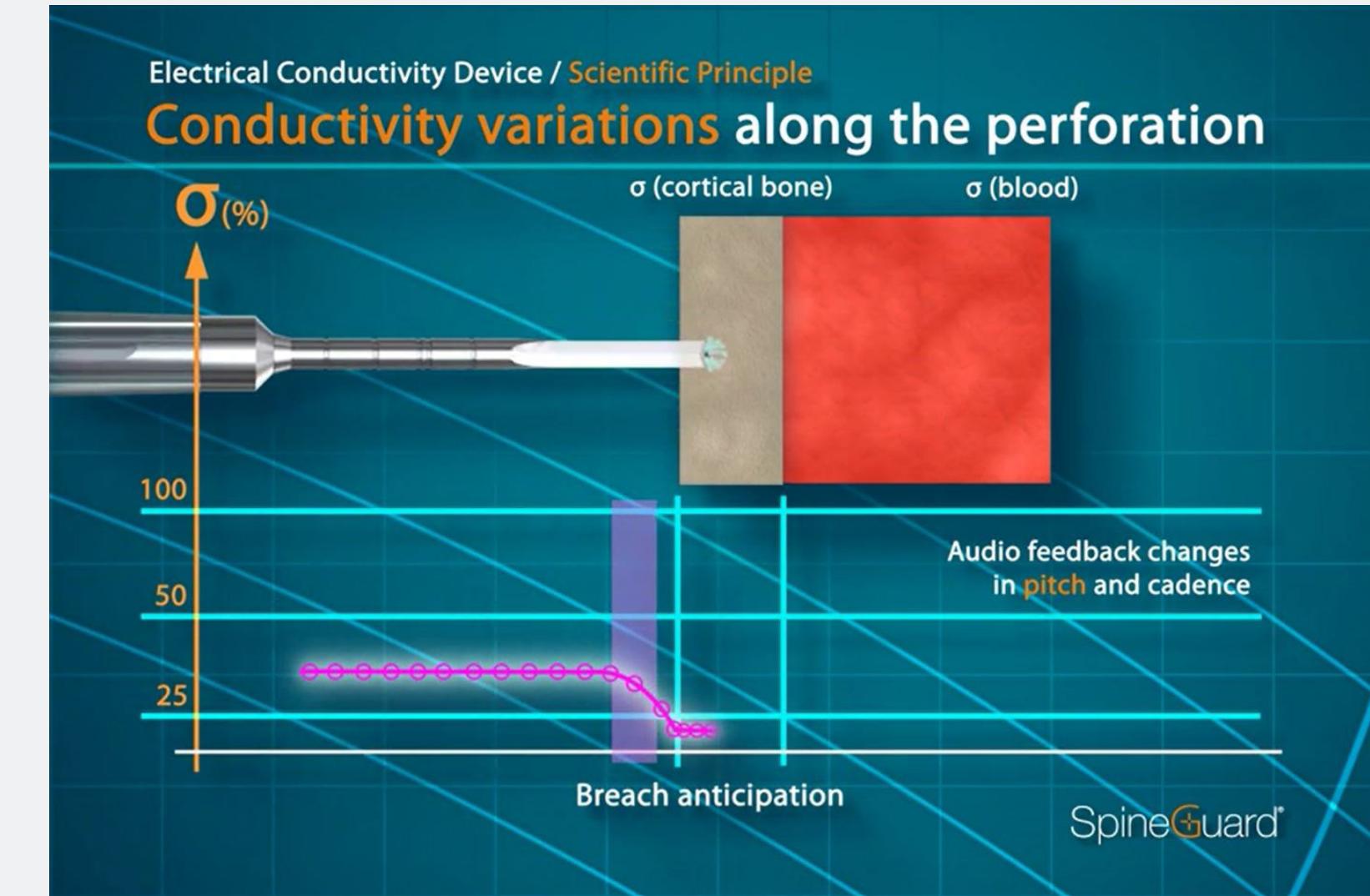Low pitch, low cadence

Tip in soft tissue/blood

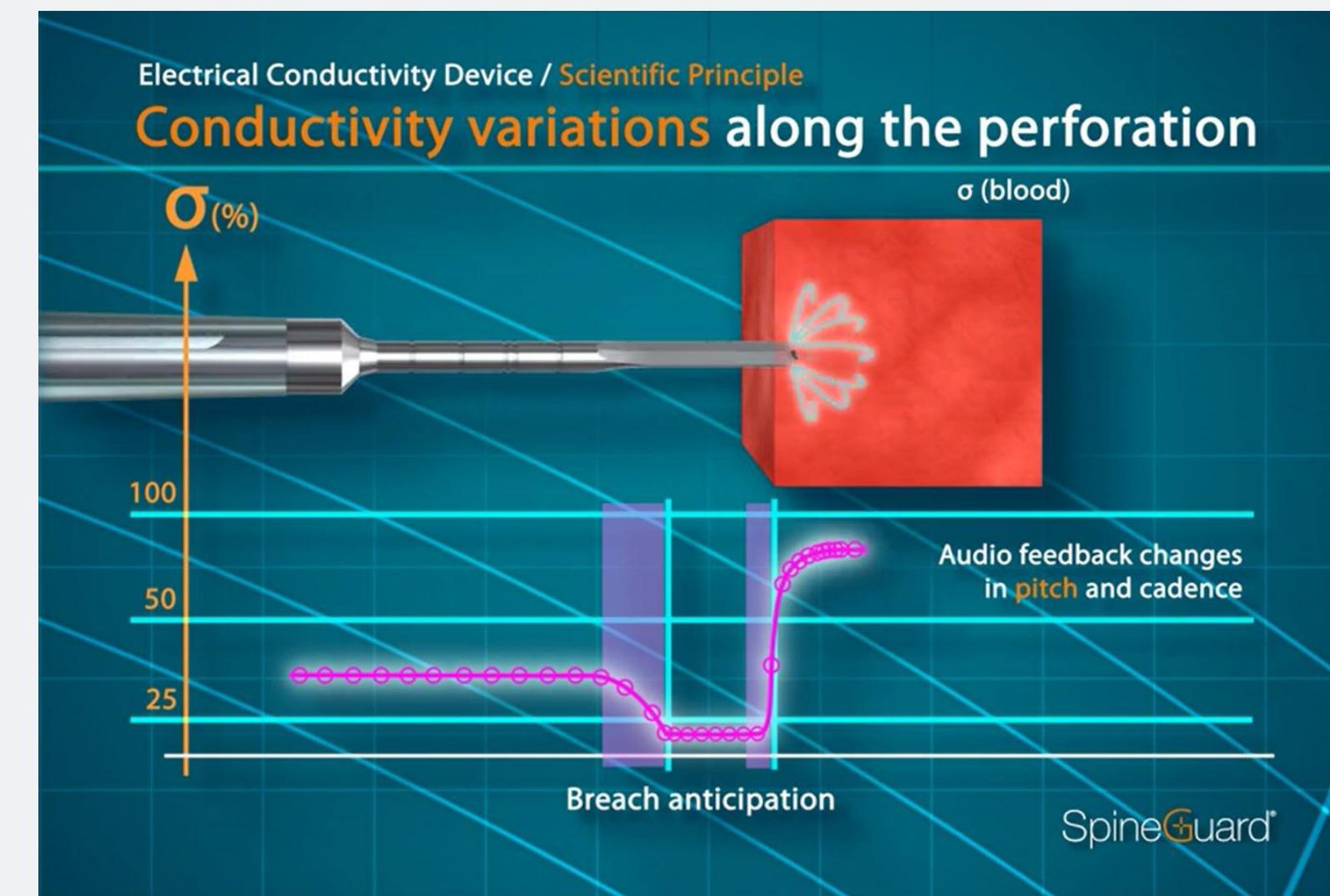High pitch, high cadence

What is DSG?


DSG Connect App objective

SpineGuard developed the DSG Connect App for surgeons to have a visualization of the audio feedback from the DSG technology during pedicle drillings.

The DSG Connect App allows collection and visualization of the signal measured at the tip of DSG Connect devices. The app allows for this data to be recorded, extracted and analyzed.

Note : The DSG Connect App is not intended for medical decision making.

Interface Global Overview

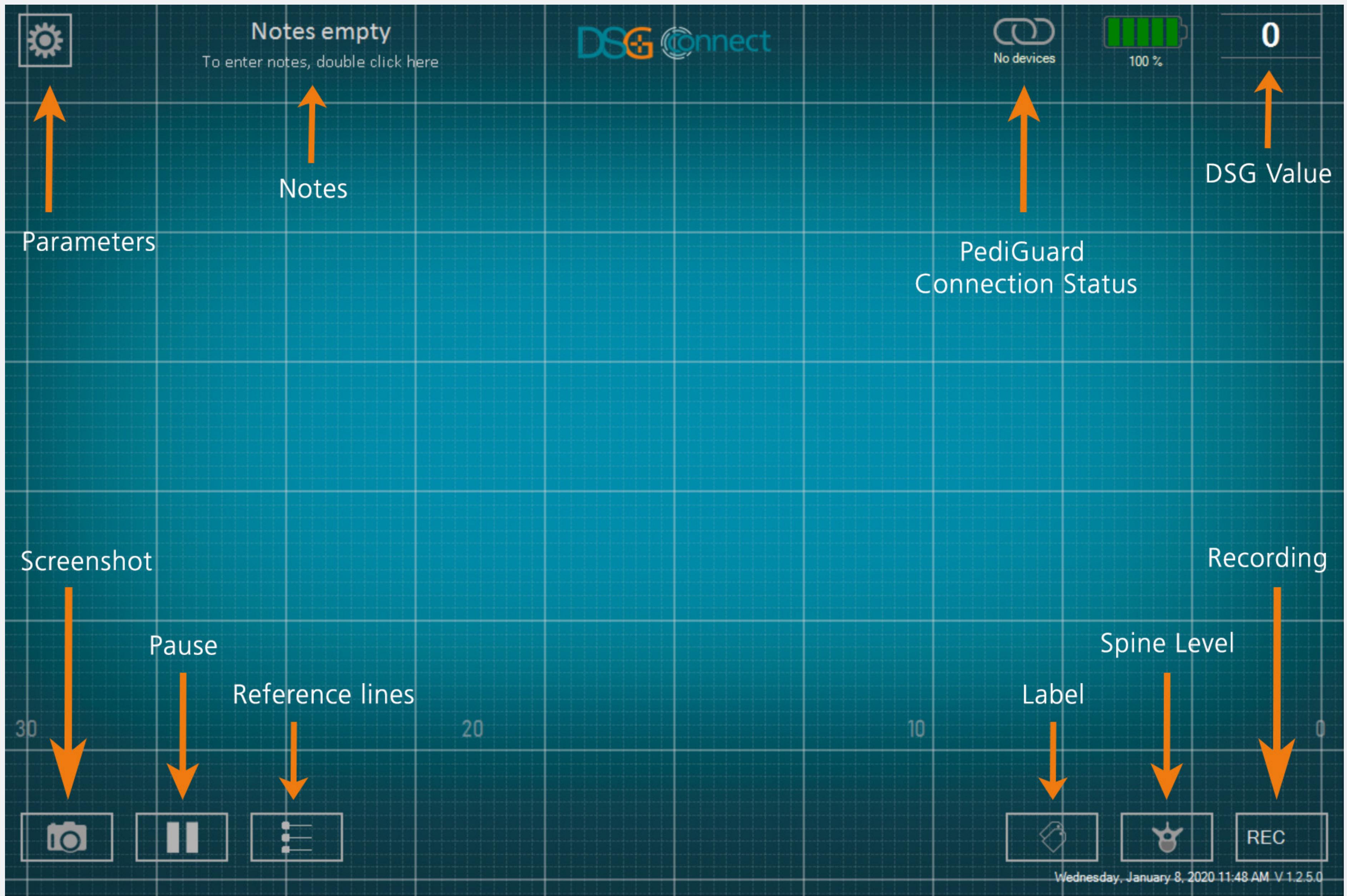
Information area

Graph area


Feature Control area

DSG Connect Device


How to connect the DSG Connect device to the tablet?

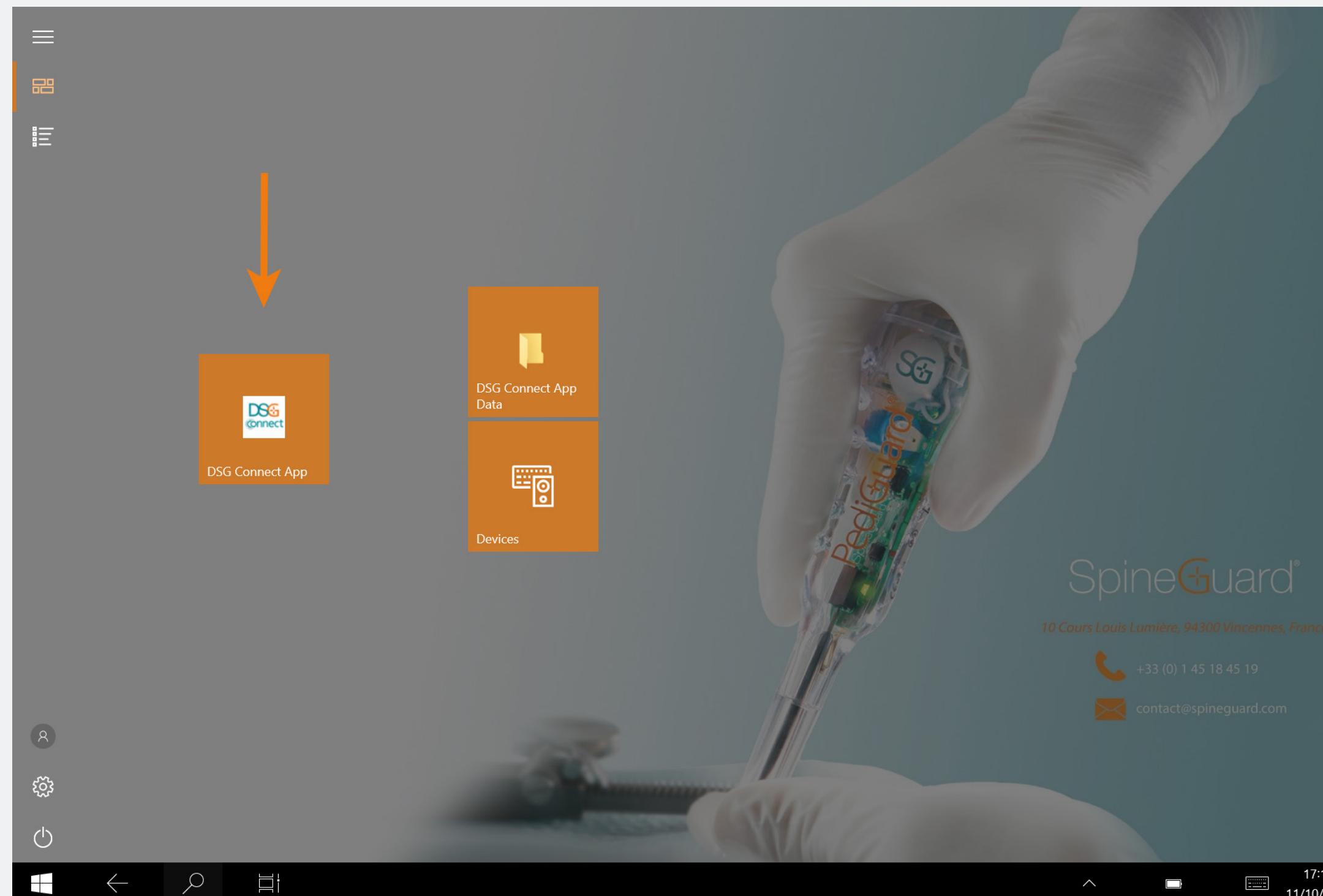- Turn the tablet on and launch the DSG App.
- Then, pull out the activation tab of the DSG Connect device. An automatic connection will be triggered between the DSG Connect device “on” and the tablet.
- The graph should start to stream and the DSG value should appear in the top right corner of the screen.

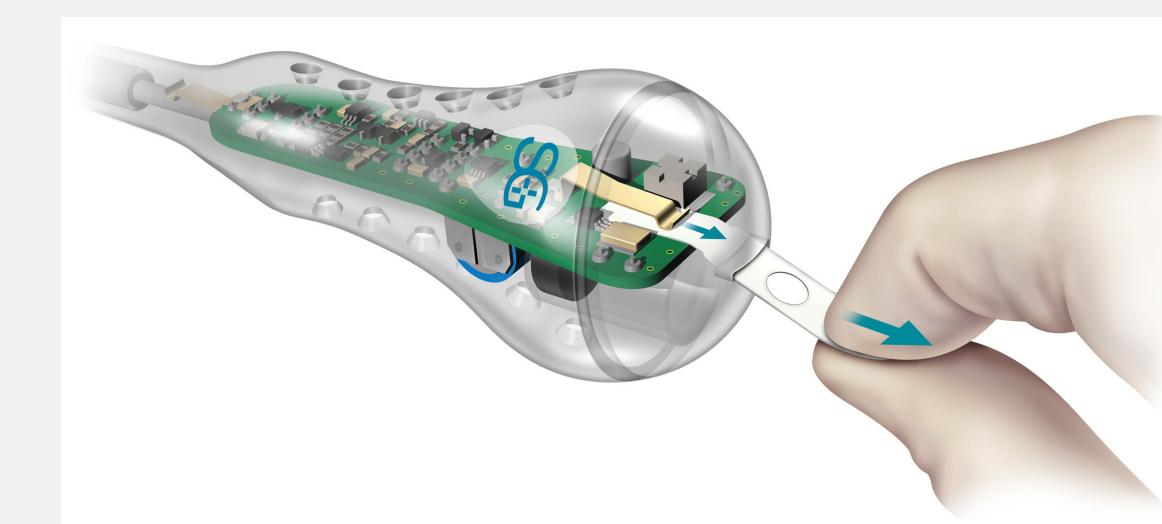What if the connection fails?

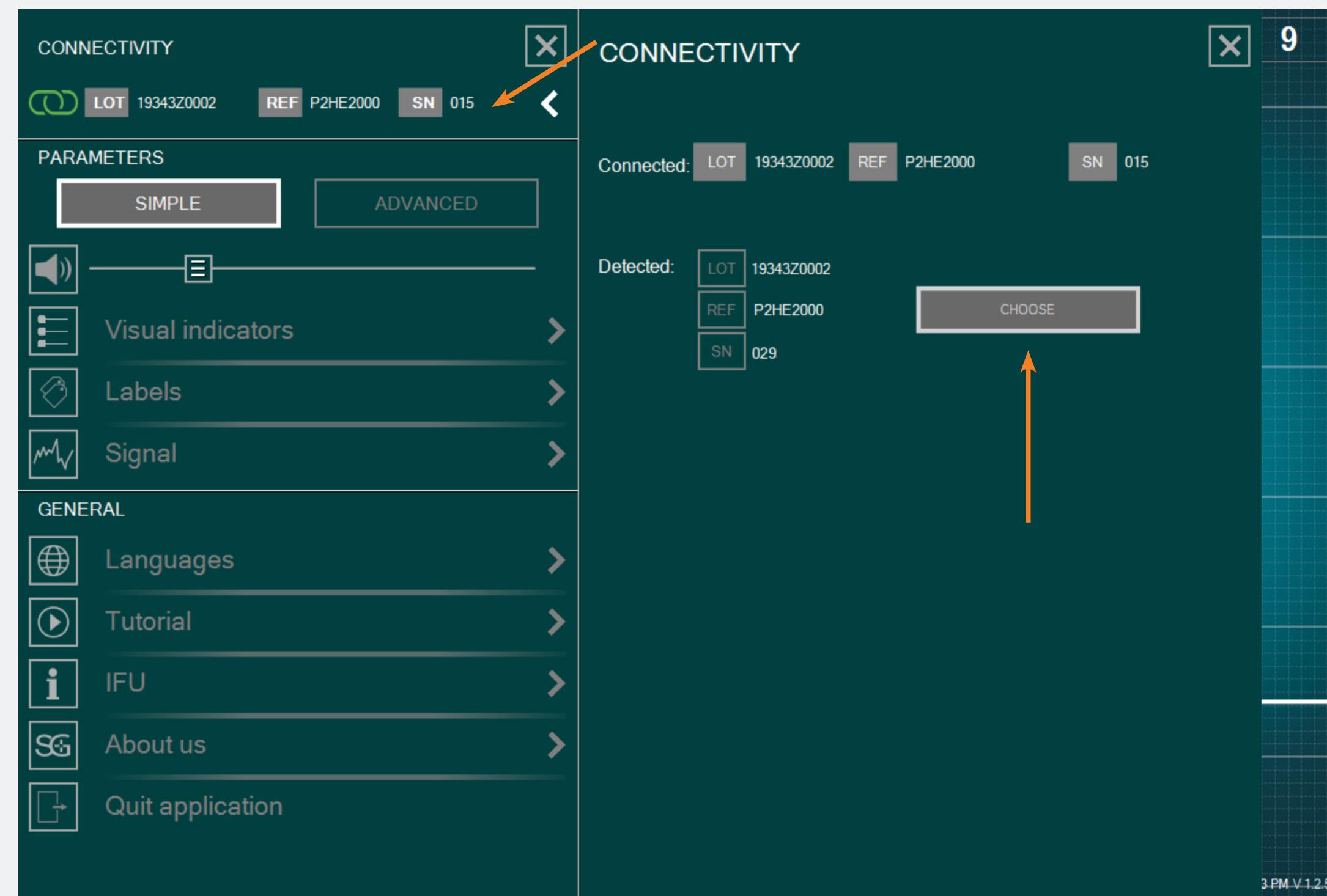- Go to the window connection in the parameters menu, and choose one of the detected DSG Connect device you would like to be connected with.

Simple Mode

The Simple Mode is a core function in which the evolution of the conductivity data along time is graphically represented by a curve on the screen. It allows you to have a visual feedback with basic features (reference lines, pause mode and screenshots). You can also personalize the Empty Field “Note”.

Signal (Graph)


What does the graph represent?

- While drilling, you will be able to visualize the data collected by the DSG Connect device.
- The conductivity data are graphically represented by the white curve on the screen.

Simple Mode

How to understand the visual signal?

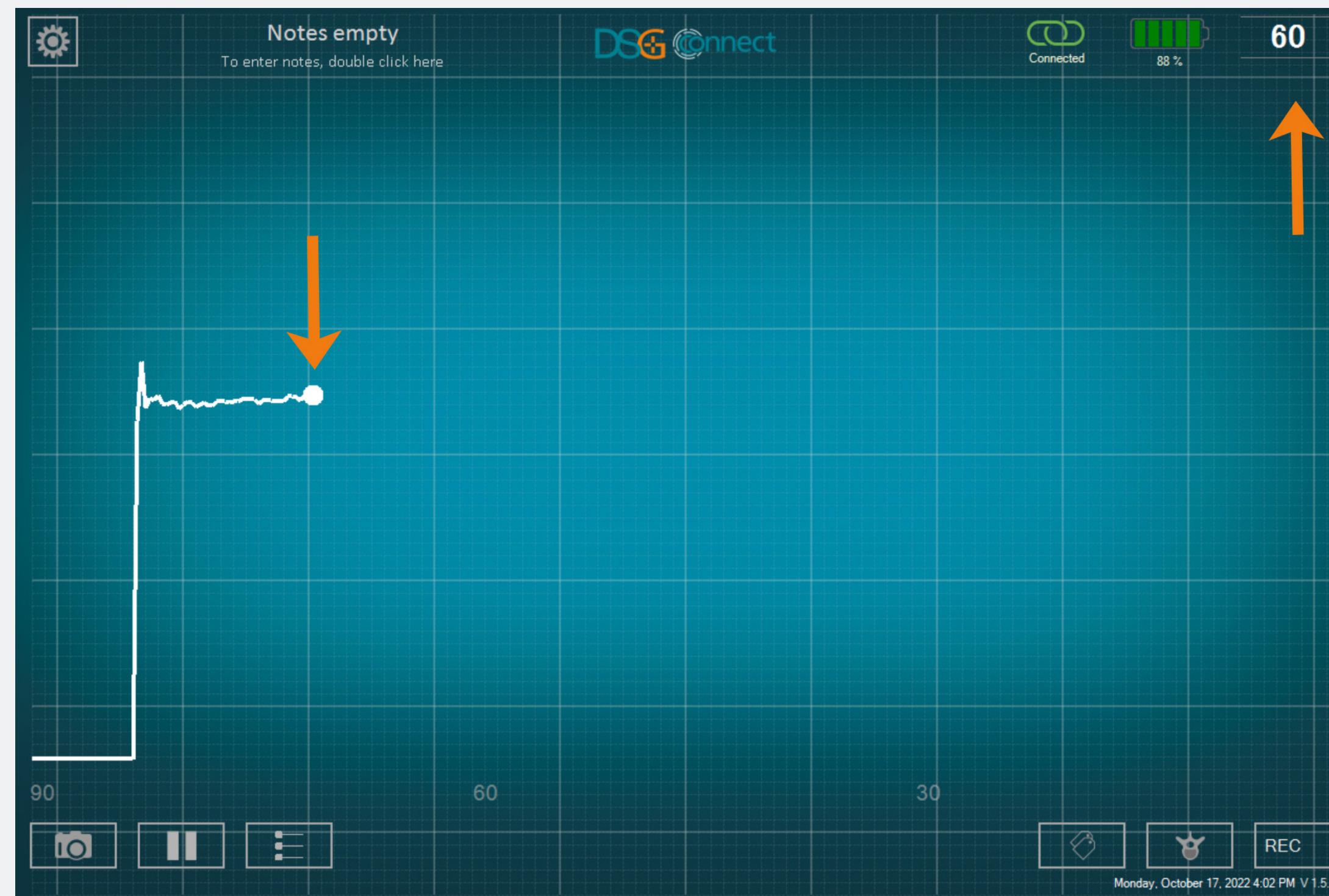- The comprehension of the signal should be made as follow: the higher the signal, the higher the conductivity and the lower the signal, the lower the conductivity.
- The white point on the graph and the number in the top right corner represent an image of the conductivity measured at the instrument's tip.
- The graph will start generating from left to right and once the full screen is filled, the entire curve will move from right to left;
- You can adjust the width of the time window in «Signal» in the parameters menu.

Simple Mode

Empty Field - Notes


What is the «Notes» feature?

- This feature allows you to enter any information deemed necessary.
- It is optional, the app can be used and data can be recorded even if the Notes field is not filled.
- If you enter a Note in the «Title» Field, the report at the end of your surgery will be recorded in a folder named accordingly.
- This text field must not contain personal data and/or health data allowing patient identification.

Simple Mode

How to add Notes?

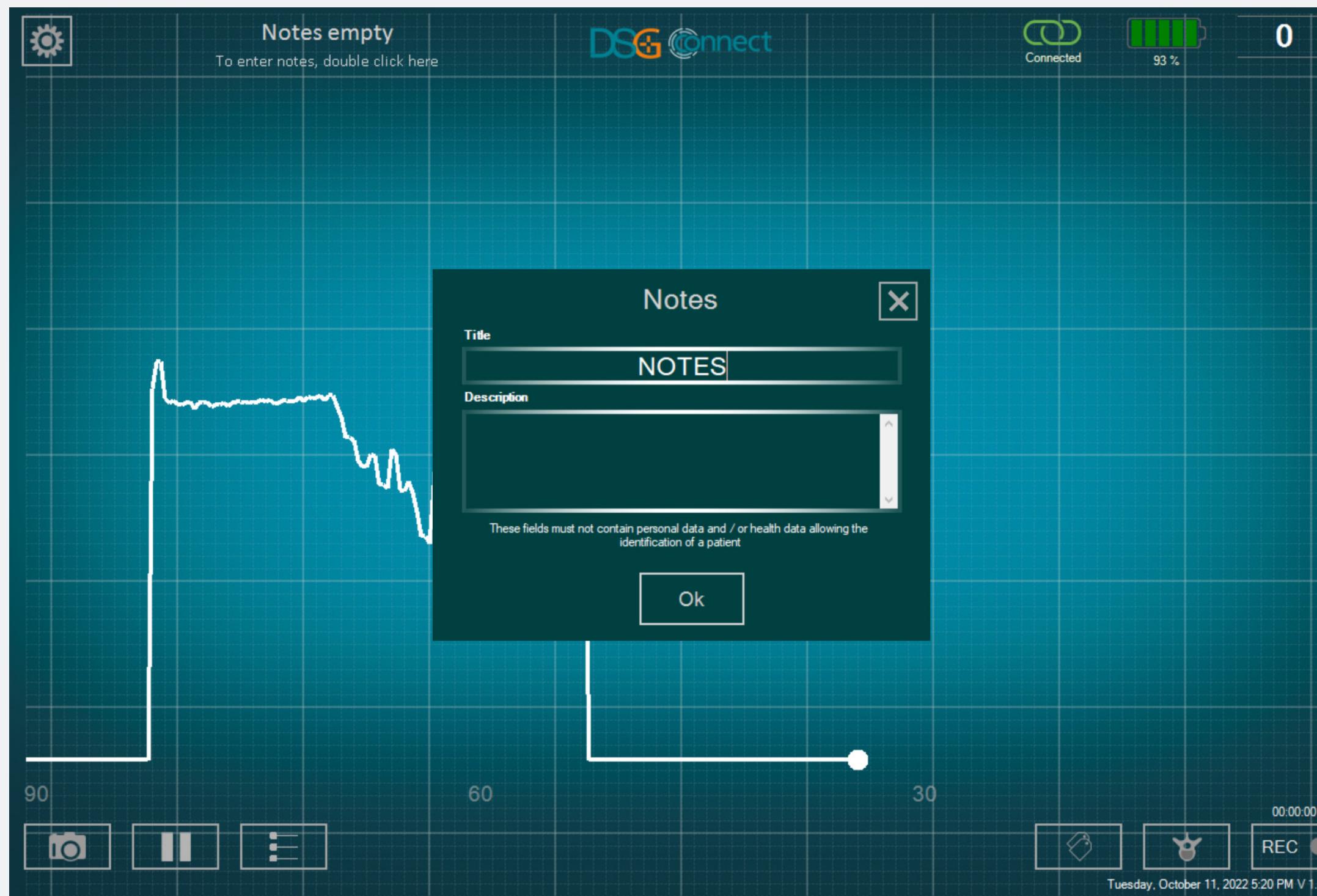- Double click on “Notes empty”.
- Fill in the corresponding empty boxes.
- The content of the note fields is maintained all along the surgery and is only reset when the app is quit. It is of course always possible to edit the content of the note during the surgery but not live during a recording.

Simple Mode

Reference lines

What is the «reference lines» feature?

The reference lines feature allows you to position a landmark on the graph which you can adjust according to your drilling.

Simple Mode

How to activate the «reference lines»?

- To activate the reference lines feature, click on the button: .
- When the feature is active and during use, the button will become bright.
- You can position reference lines by dragging the desired color cursors up and down.
- To add reference lines, drag one of the buttons on the left panel  (on top or bottom), and position the reference lines where you want on the scale.

Simple Mode


How to adjust «reference lines» position?

- Drag one of the buttons  at the desired location.
- The value of the reference line is viewable on the right side of the screen.

Simple Mode

How to delete «reference lines»?

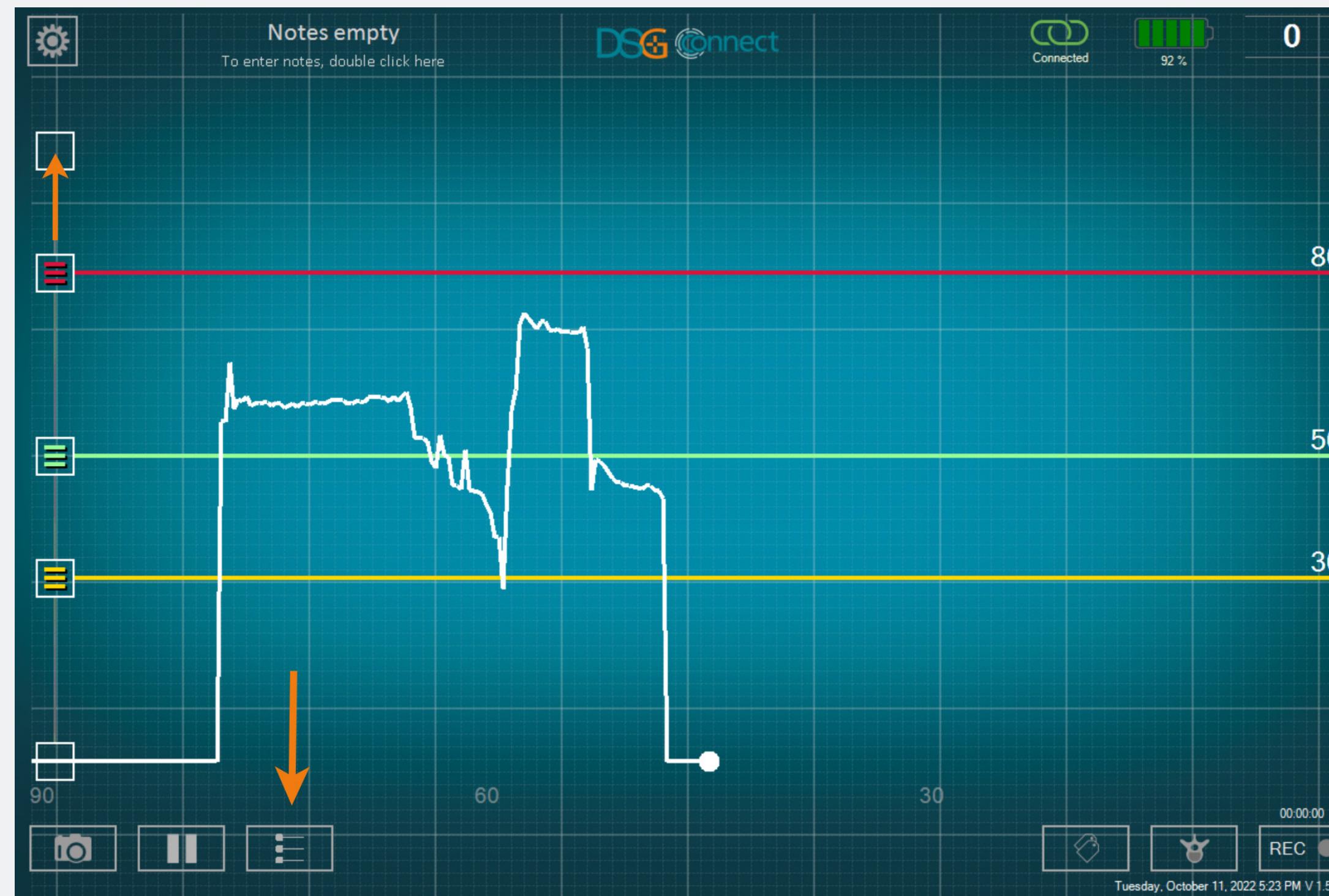- You have the possibility to remove the reference lines individually by sliding the corresponding button up or down the scale until the end of the selection zone.
- If you want to hide all the reference lines, click on the button .

Screenshots


How to use the «Screenshot» feature?

- You can take screenshots at any moment of the surgery.
- Screenshot your graph if needed by clicking on the screenshot button

Simple Mode

Where can I find the screenshots I took?

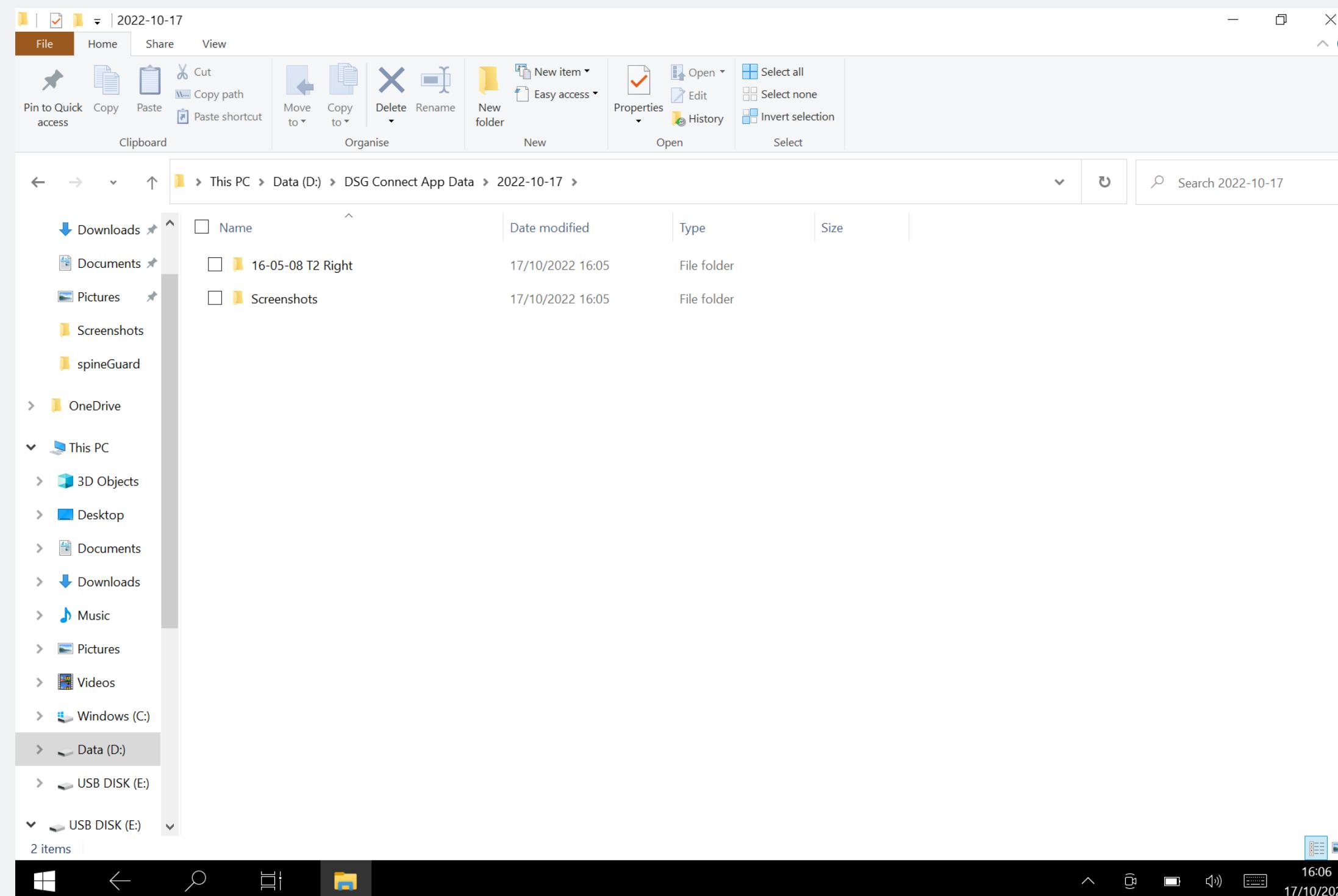You can see your screenshots at the end of the surgery in the DSG Connect App data. See Export-Report section for more details.

Simple Mode

Pause

How to use the «Pause» feature?

- Click on the Pause button: .
- The screen will freeze and a notification will appear.


To go back to normal, click once again on the Pause button: .

Advanced Mode

The Advanced Mode is mainly used for data recording purposes. It has additional features compared to the Simple Mode such as Spine Level Selection, Recording, Labels.

Pedicle Selection


What is the “Pedicle Selection” feature?

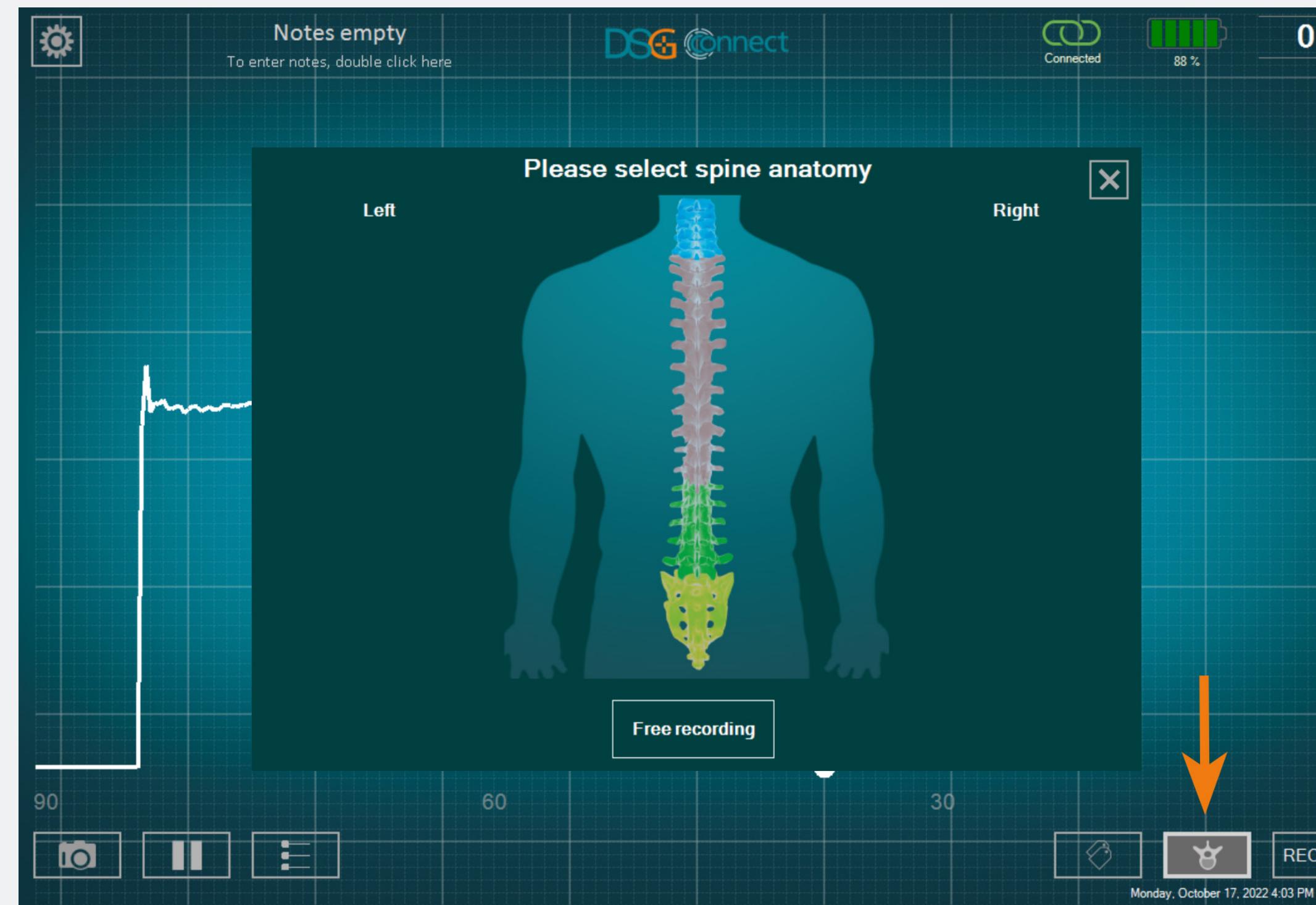- The Pedicle Selection feature allows identification of the pedicle for each drilling record.
- You have to select your pedicles to be able to use the record mode.
- The data are recorded by pedicle drillings to facilitate the review at the end of the surgery.

Advanced Mode


How to use it?

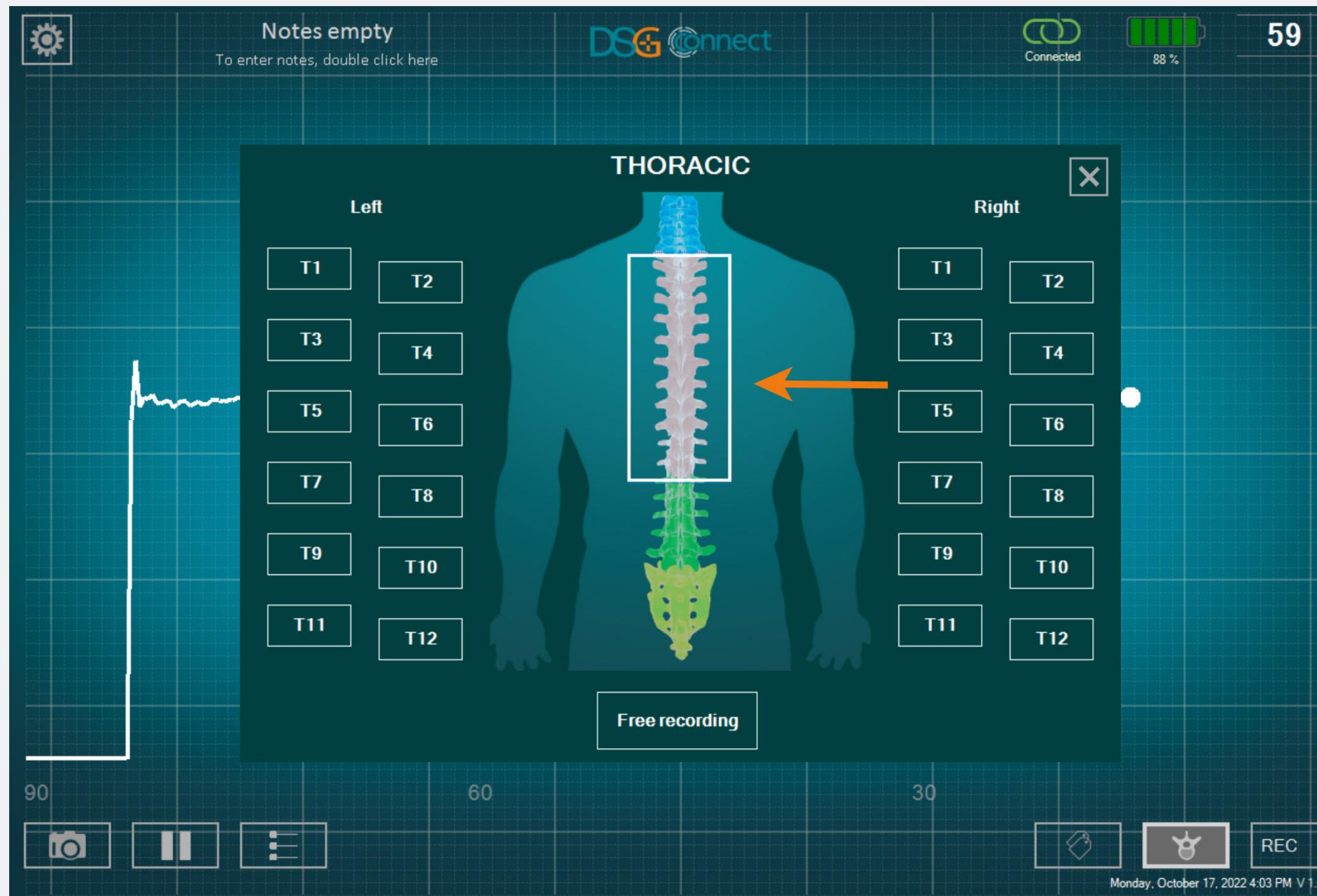- Click on the Spine Level button: . A window with a spine appears.
- Select the level of the spine you want to drill.
- Then, select the corresponding pedicle. The state of the button will change from to .
- After having selected your pedicle, the window automatically closes.
- The pedicle you selected is indicated on top of the Spine Selection button.

Advanced Mode


What if I want to select a pedicle I have already selected, and already recorded the drilling?

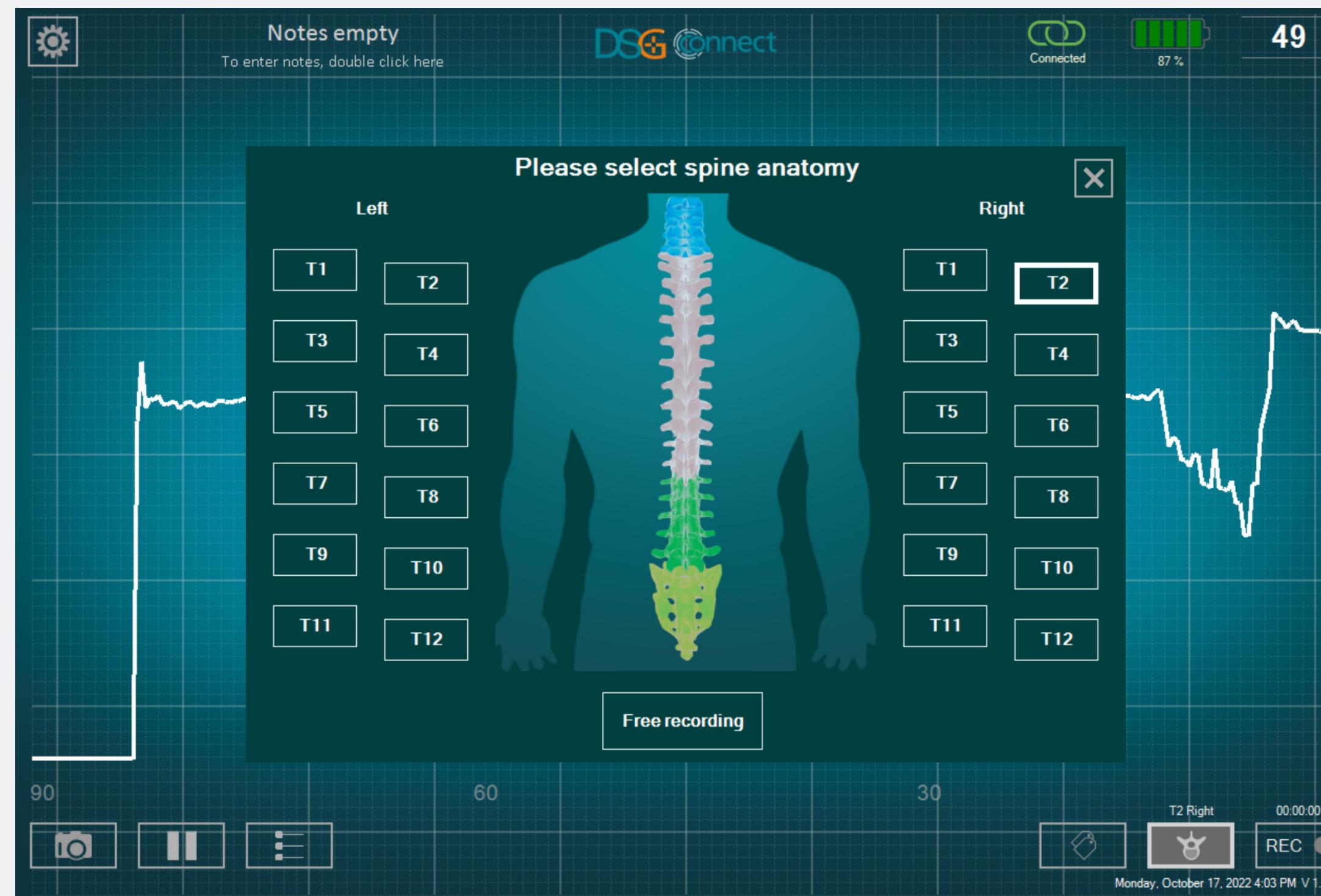Click a second time on the corresponding button. The button (as it has already been selected) will be shaded: **T2**.

Advanced Mode

What is “Free Recording”?

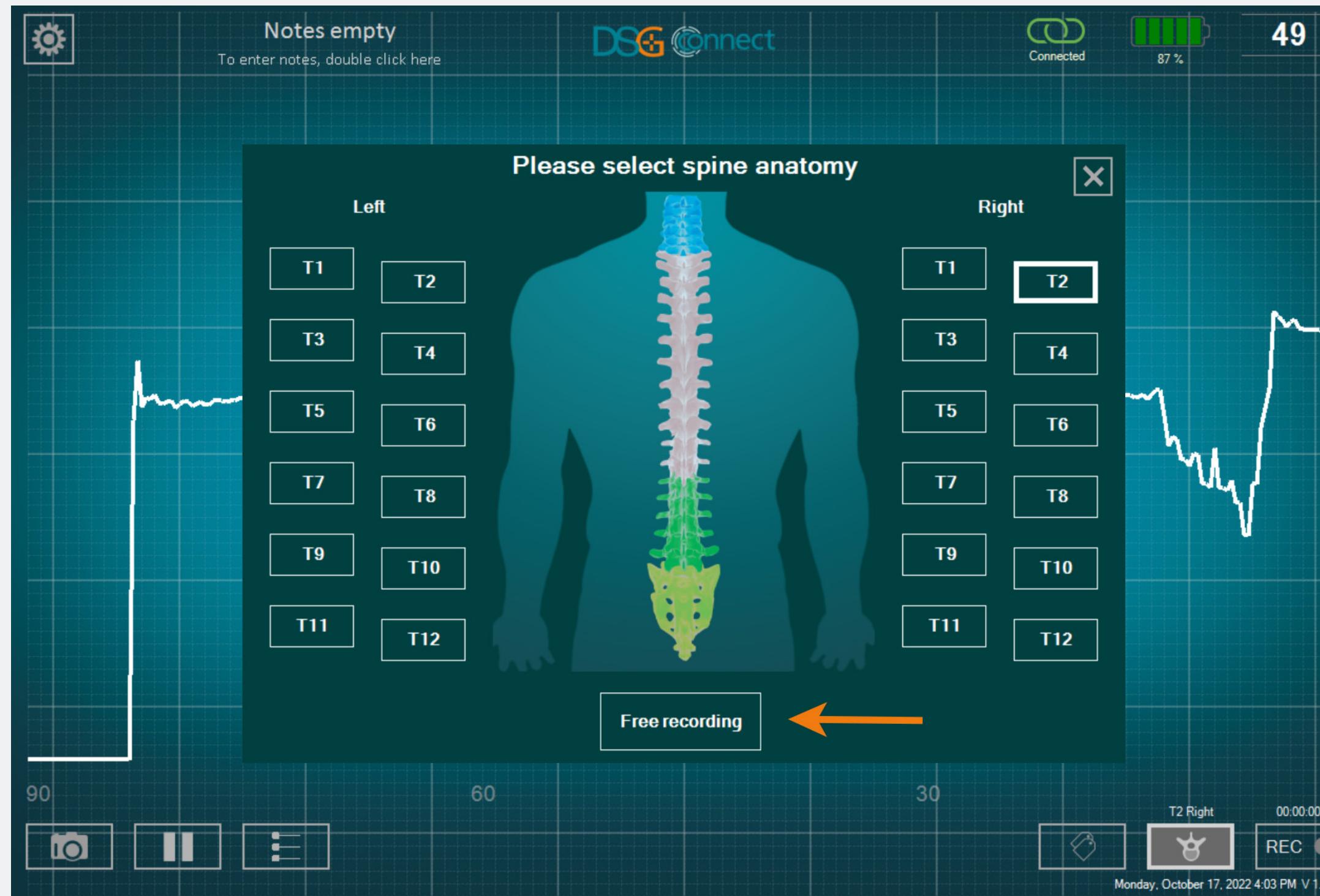
If you do not want to record the data per pedicle, you can select “Free Recording”. It works the same way as the pedicle selection.

Note: Recording is only possible if a spine level or free recording is selected.

Recording


How to use it?

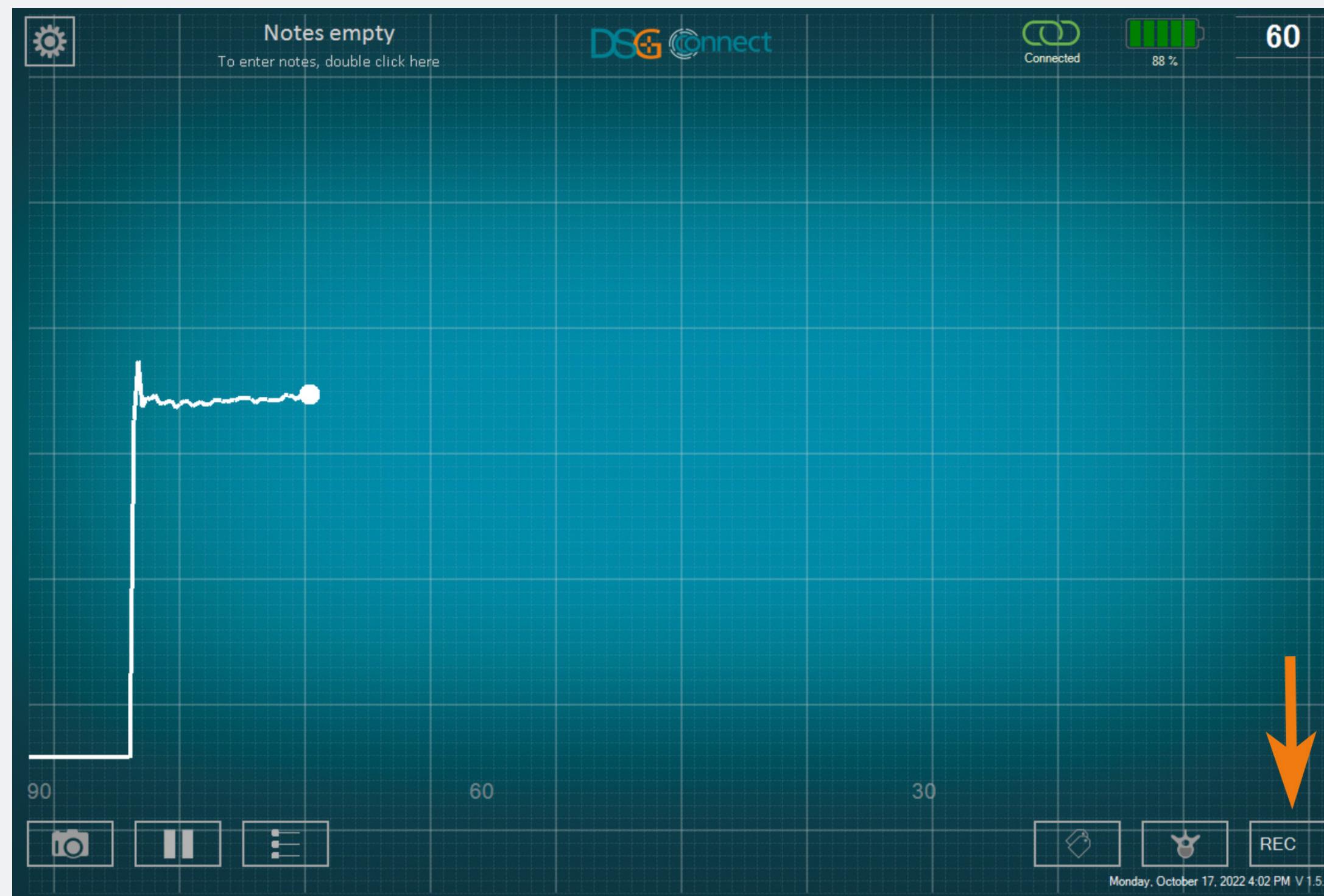- After having selected your pedicle or Free recording, a grey circle in the REC button will appear: . It means you can start to record your drilling.
- Click on that button to record. The grey circle will turn red and will blink, it means that you are recording: . The time-indicator will inform you of the duration of the recording.
- At the end of the drilling, If you want to stop your recording, click once again on the button.

Note: When a recording is started, the graph restarts at 0.

Where can I find my recordings?

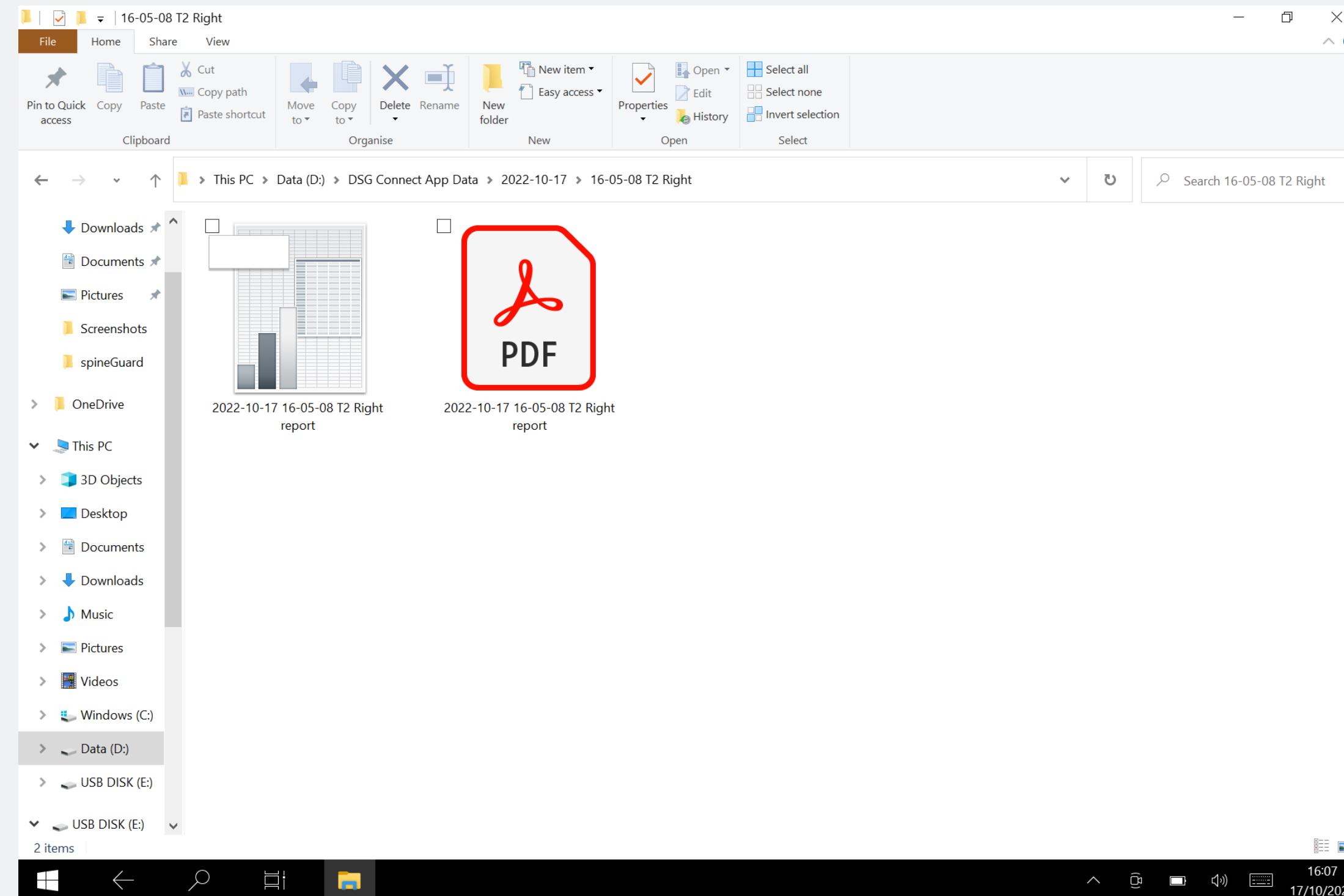You can see your recordings at the end of the surgery in a dedicated folder. One spreadsheet file and one PDF file are generated for each recording done.

Advanced Mode

Label


What is the “Label” feature?

It allows you to identify a specific event during the surgery by putting a tag on the screen.

Advanced Mode

How to use it?

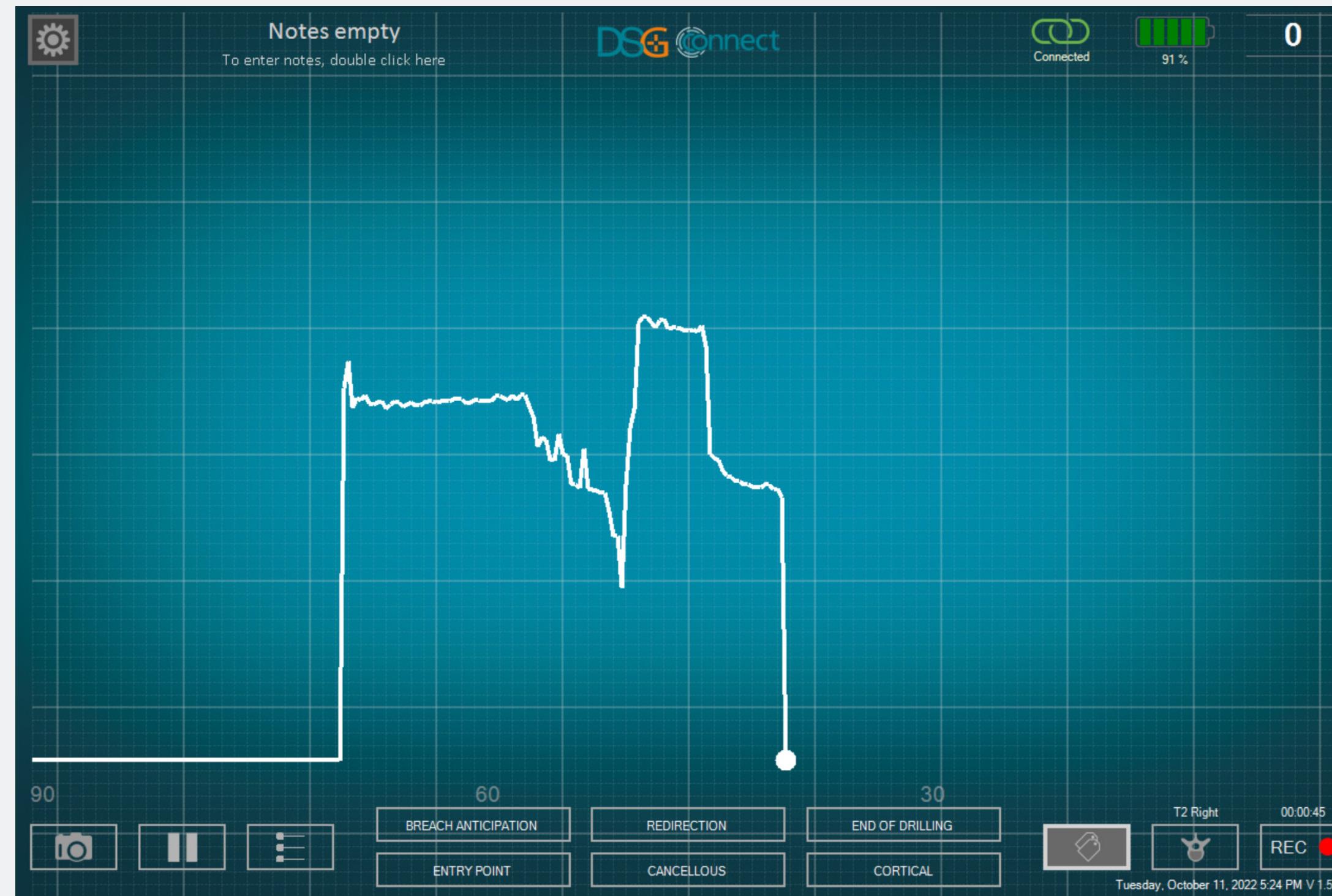- When you click on the «Label» button: , the predefined label buttons appear.
- Click on the desired button to automatically position a label on the current plot point.
- If you want to insert another label, repeat the process.

How to reposition a «label»?


Once the label has been placed on the graph, it is possible to reposition it.

- Click on the concerned label. When the label is selected, it becomes bright.
- Click on the desired location on the graph. The label will be automatically repositioned.


Export-Report (data review)

Where to find my screenshots and my recordings?

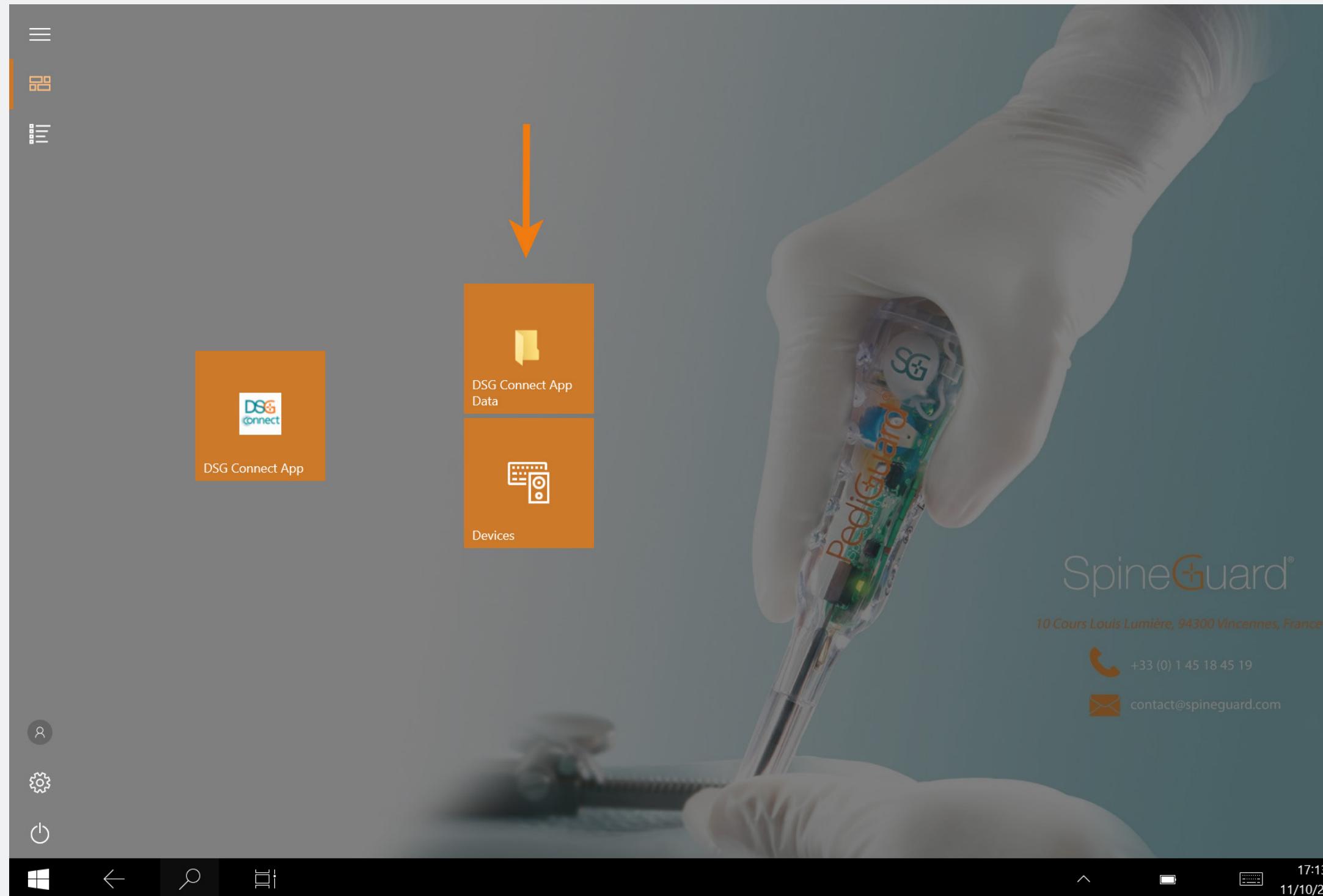- Your screenshots and recordings are available after the surgery and are available by clicking on the DSG Connect App Data tile shortcut on the homepage. They are stored in separate folders that are classified by date as follow:
 - Name: YYYY-MM-DD
 - Example: 2020-01-07
- If you entered a Note, the folder will be named as such:
 - Name: YYYY-MM-DD Note
 - Example: 2020-01-07 NOTES 1

Advanced Mode

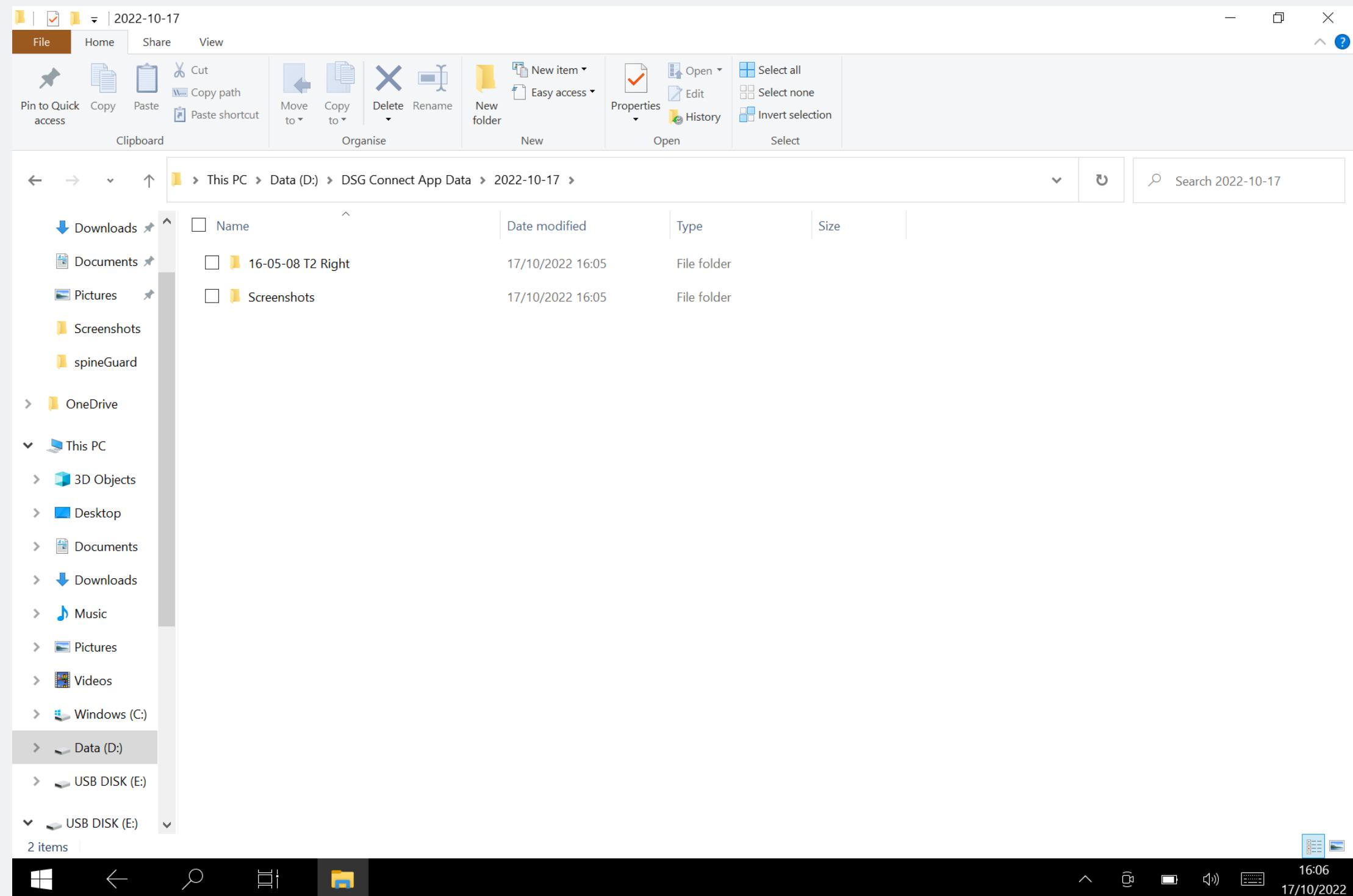- Recording:
 - Format: CSV and PDF
 - Name: YYYY-MM-DD HH-MM-SS Selected_Pedicle report
 - Example: 2020-01-07 11-44-50 T2Right report


- Screenshot taken with no recording ongoing:

- Format: JPG
- Name: YYYY-MM-DD HH-MM-SS
- Example: 2020-01-07 11-43-24
- Location: folder of the day

- Screenshot taking during the recording:

- Format: JPG
- Name: YYYY-MM-DD HH-MM-SS Selected_Pedicle screenshot
- Example: 2020-01-07 11-44-50 T2 Right screenshot
- Location: folder of the recording


What does the report contain?

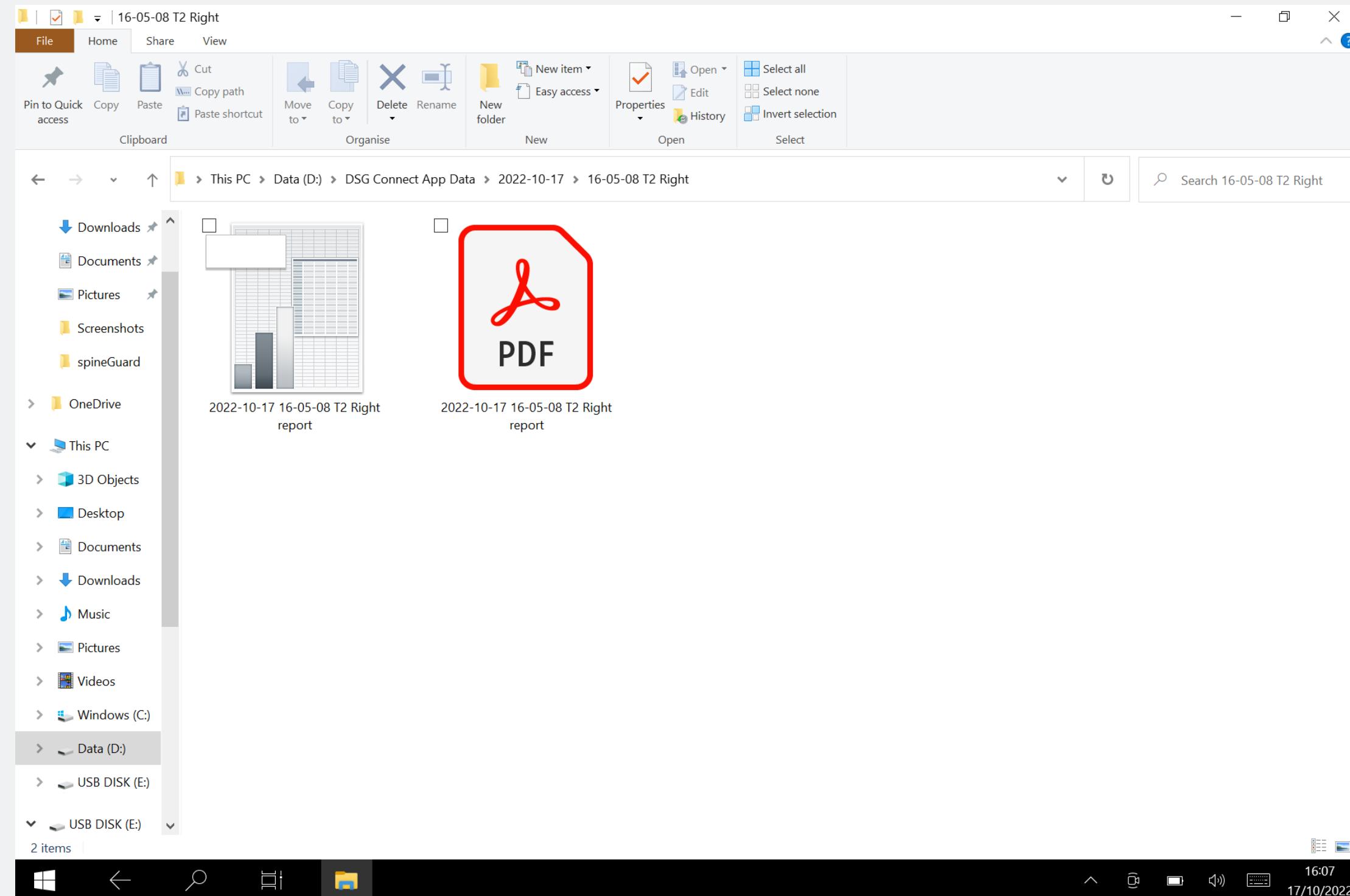
- The PDF file contains the basic data such as:
 - The device data
 - The start and end time of the recording
 - The notes added during the surgery
 - The graph
 - The labels
 - The reference lines.
- The .CSV file contains the same information but in a raw spreadsheet format.

DSG Device Connection

How to know if the DSG Connect device I want to use is connected?

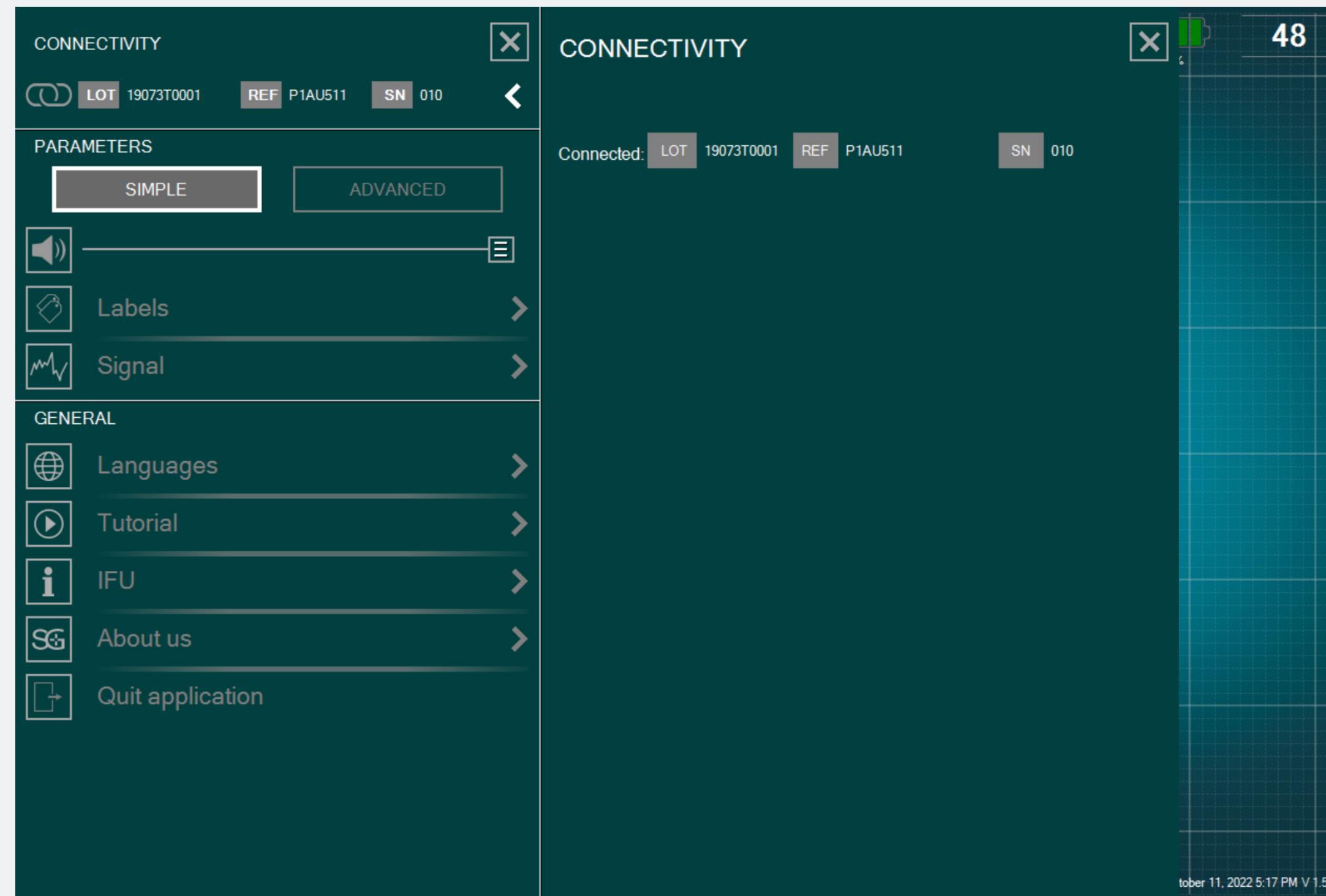
- The green connected icon means one DSG Connect device is connected to the tablet:
- The traceability information of the device (reference, lot number and serial number) is indicated in the connection window.
- If there is no DSG Connect device Connected:

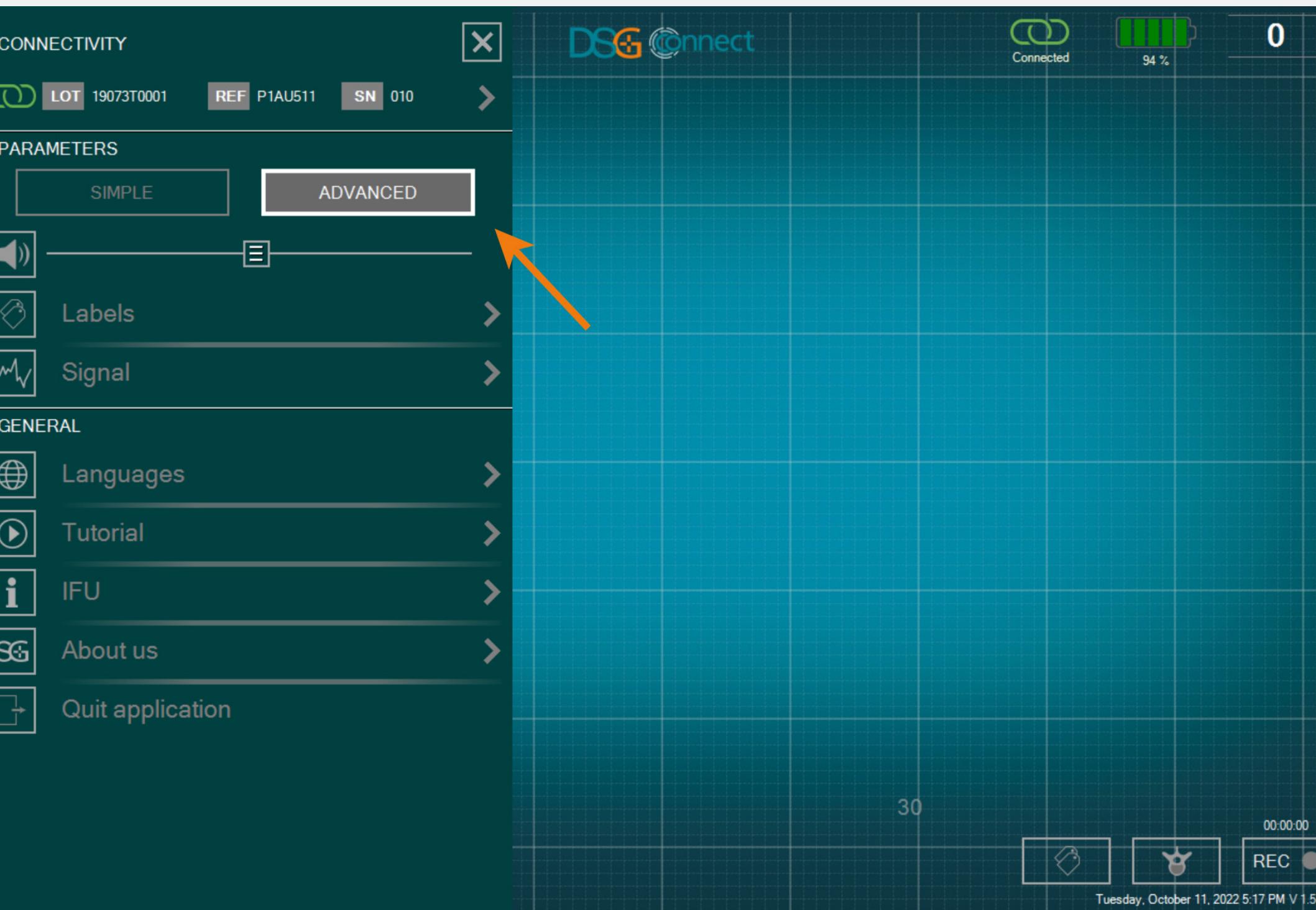
You can select either the Simple mode or Advanced mode, depending on your needs during surgery.

Parameters

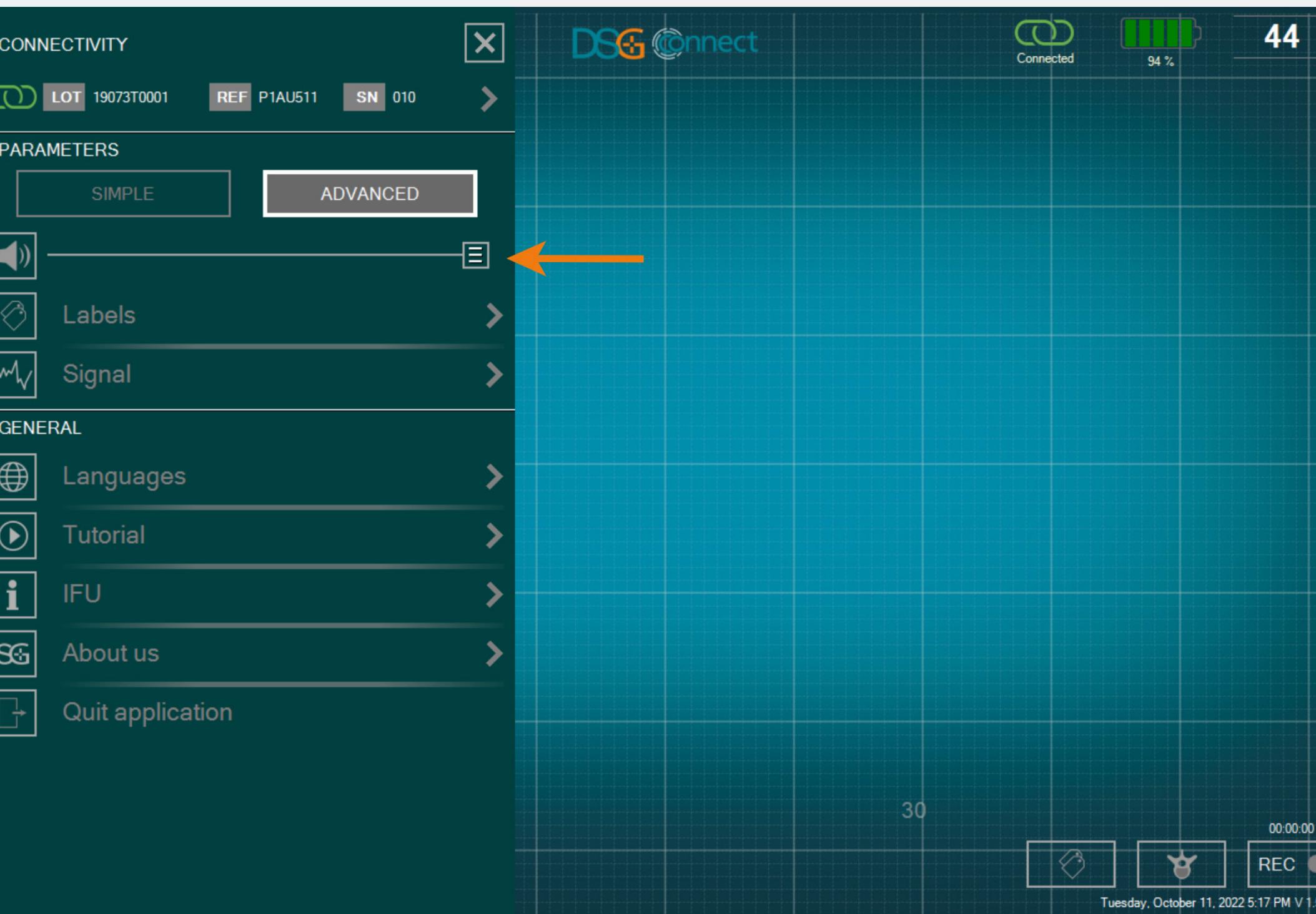
- The sound notification volume can be modified by sliding the corresponding bar from left to right.
- The sound notification will indicate:
 - Low battery,
 - Connection loss,
 - Full Memory,
 - Pause,
 - Screenshot.

Parameters

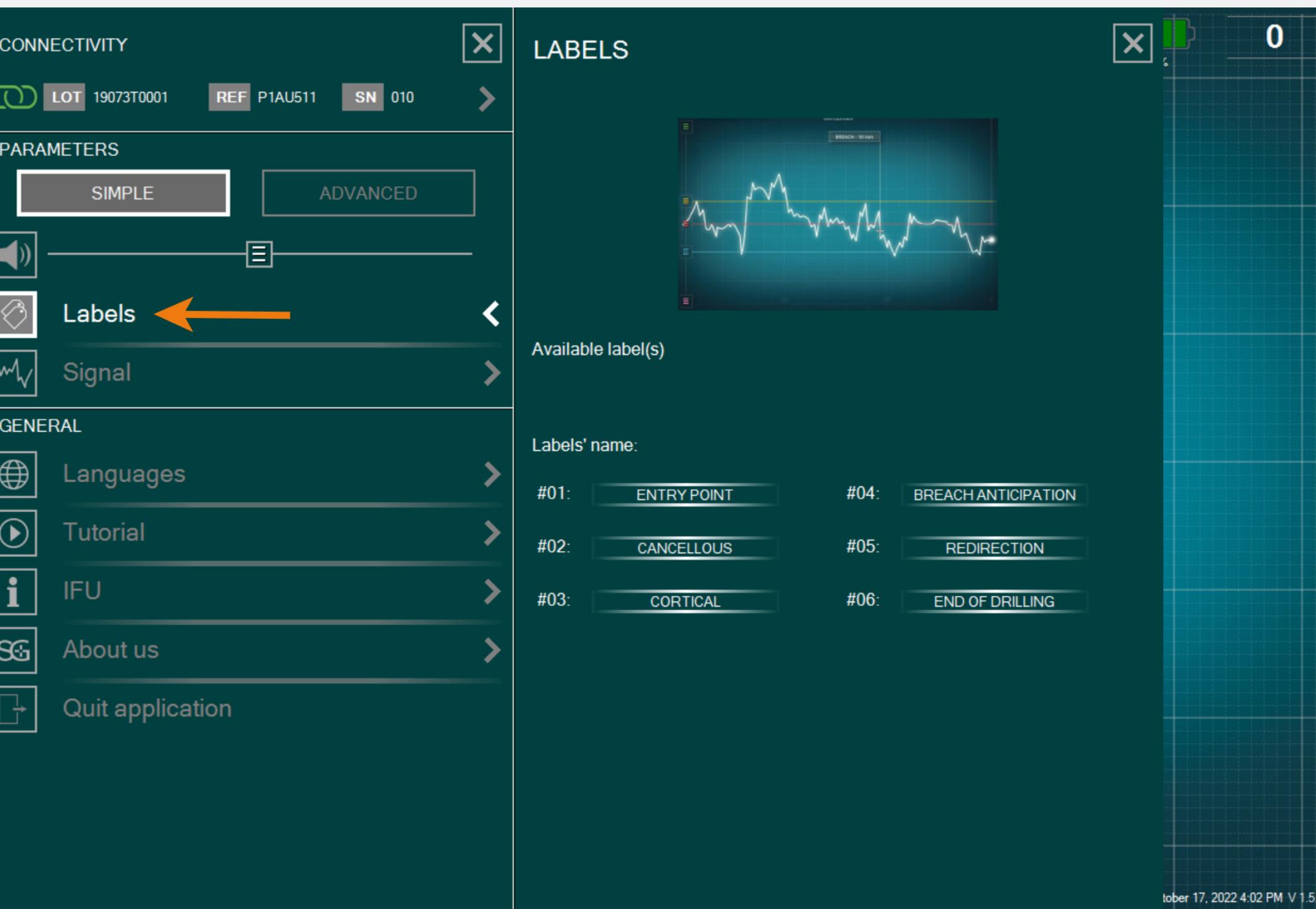
There are several pre-registered labels with the most frequently used labels by default. You can edit these labels by clicking on the text box. A keyboard will appear and you can enter your custom label (up to 20 characters).

Parameters

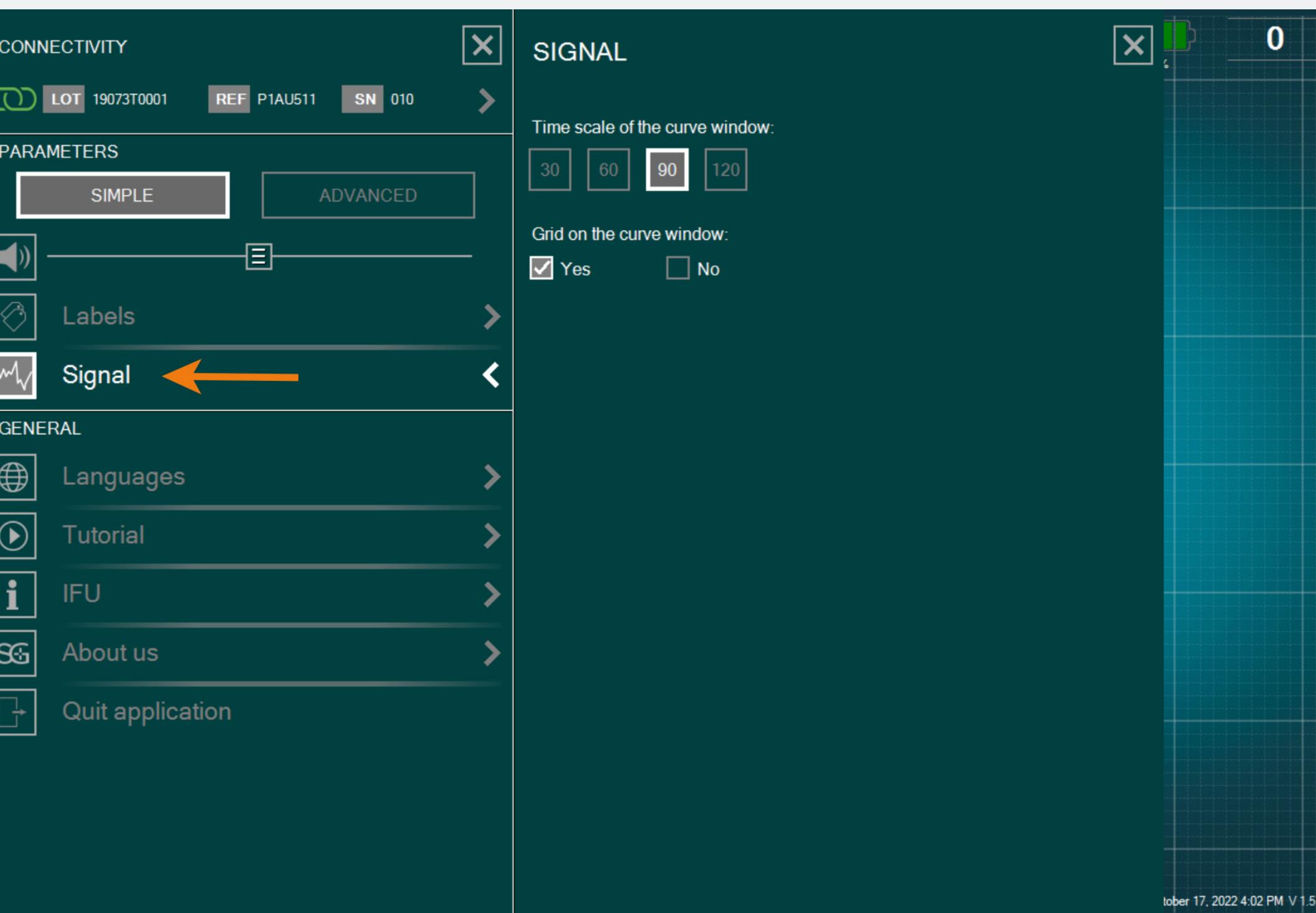
- You can choose the width of the time window to be 30s, 60s, 90s (by default) or 120s.
- If you do not want to have the grid background on the main screen, select “No” for the option “grid on the curve window”.

Parameters

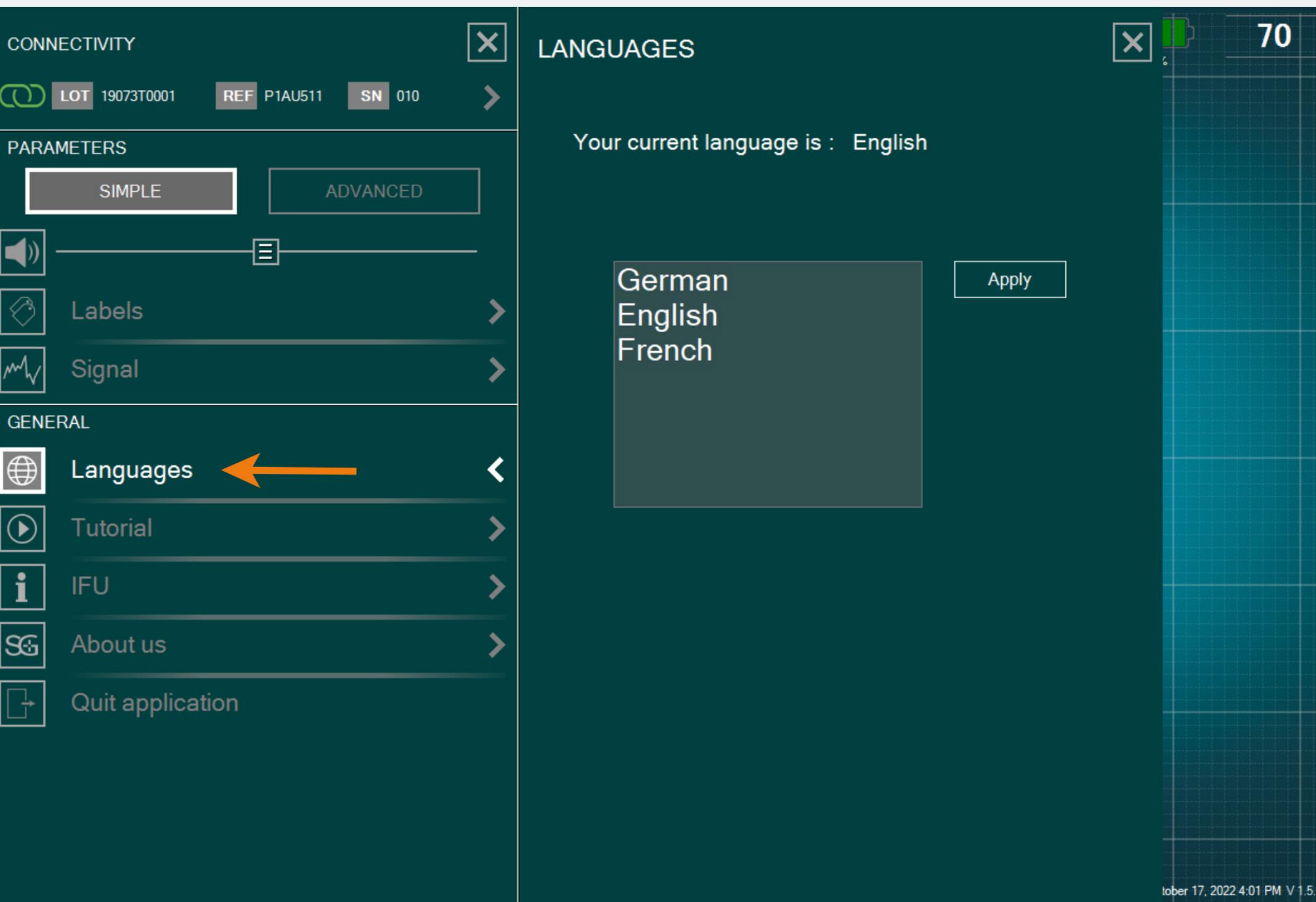
Select your language.

Parameters

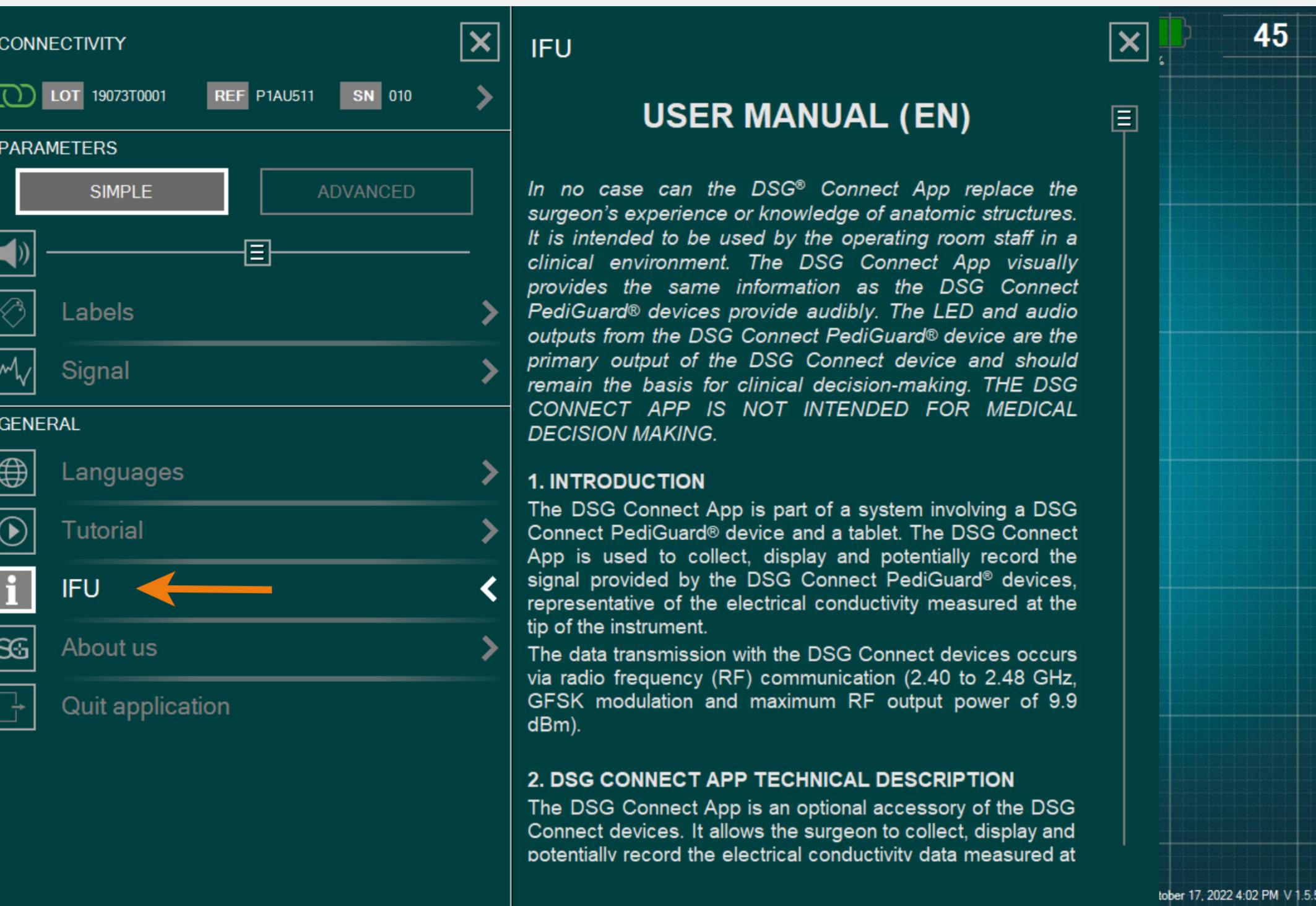
The IFU of the DSG Connect App is available in the parameters.

Parameters

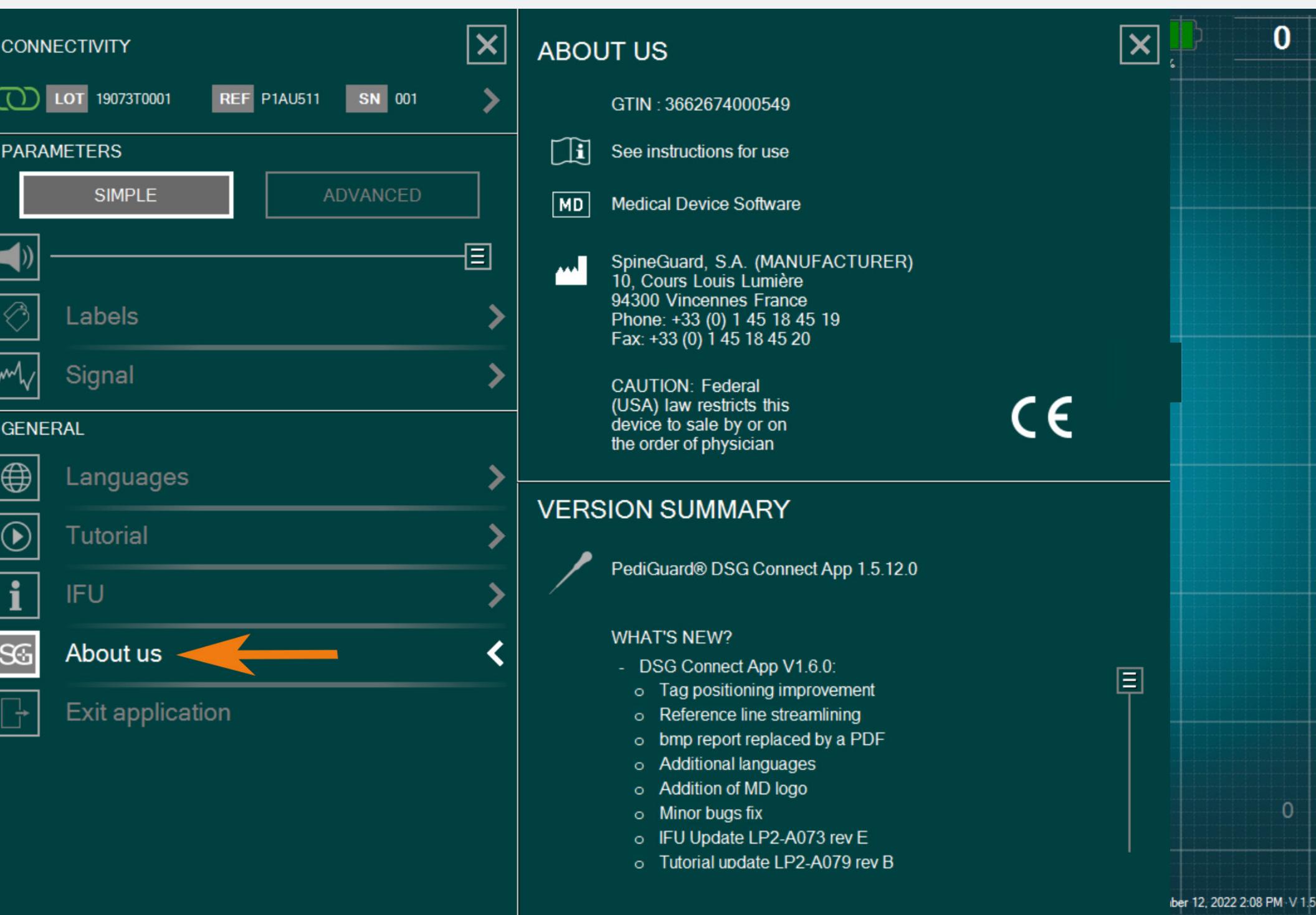- “About Us” refers to the general company information and the Version Summary.
- The Version Summary refers to the version of the DSG Connect App and current version.

Parameters

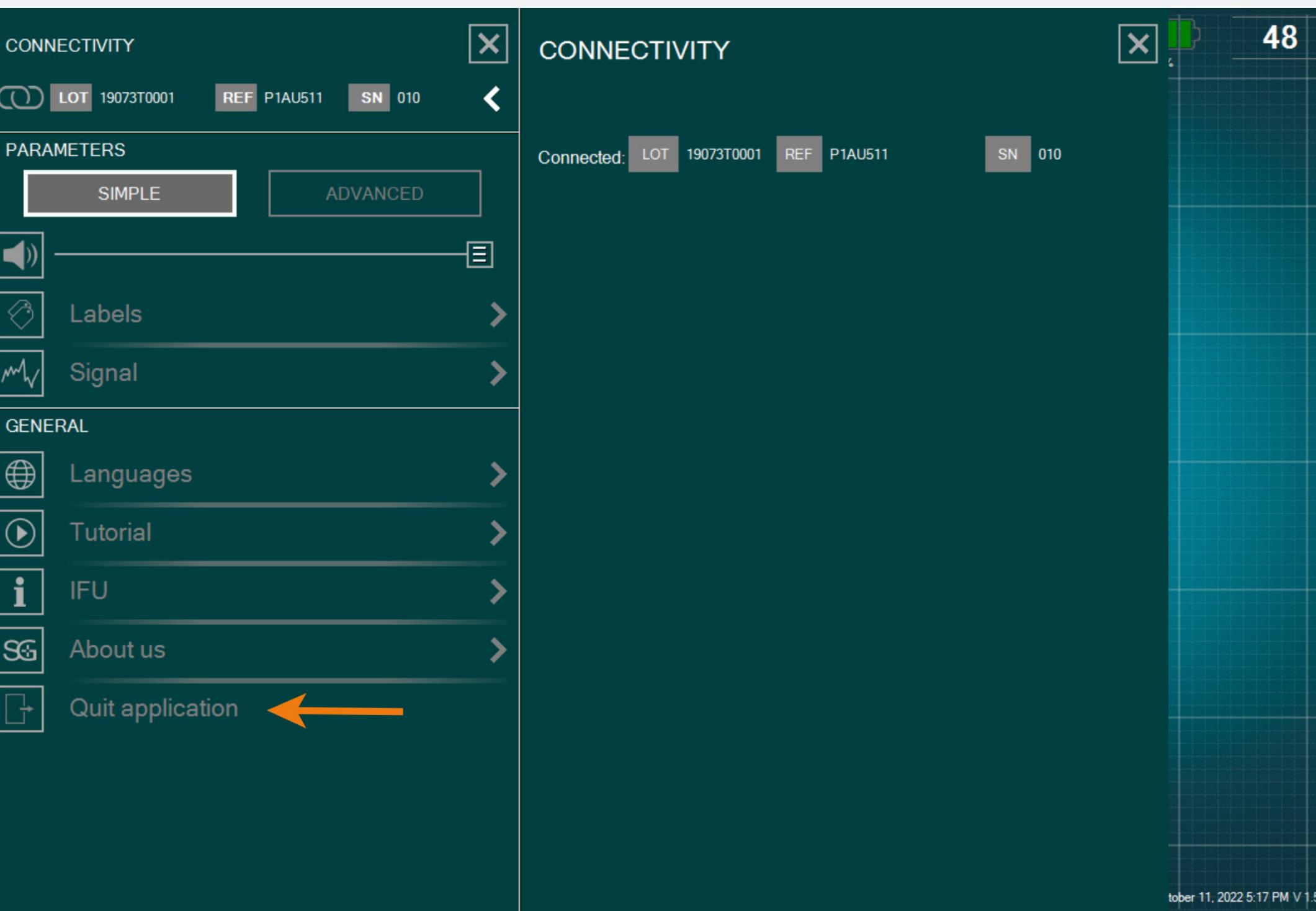To exit the App, click on the button «Quit application»:  and confirm your choice by clicking «yes».


System Notifications

- From 20%, the “Low Battery” notification will pop up.
- At this stage, we advise you to plug the power cord to the tablet to recharge the battery.

System Notifications

Connection

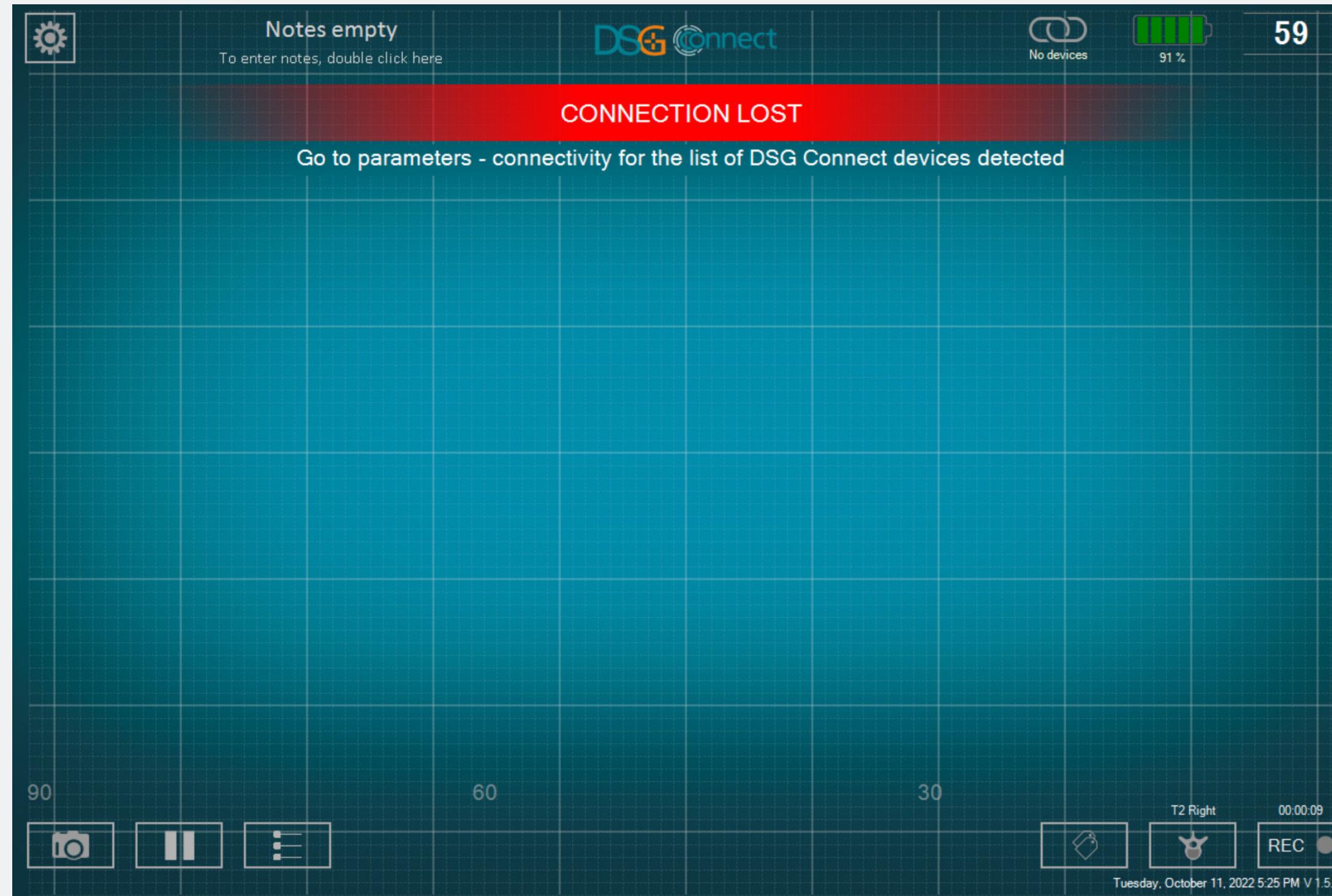
If the wireless communication is lost, go to the settings menu and select an active DSG Connect device.

System Notifications

- This notification will pop up if you have low memory on your tablet.
- It is recommended to clean the tablet memory regularly to prevent the tablet to run out of memory.
- This could result in the software slowing down, or cause it to freeze, or resulting in the inability to review live data or record additional data.

For any question, please contact SpineGuard
at contact@spineguard.com

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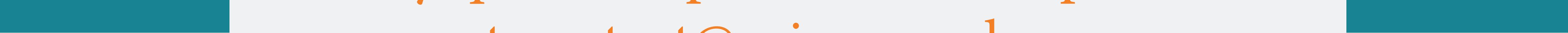Paramètres

The screenshot shows the application's settings menu. At the top left, there is a 'CONNECTIVITÉ' section with icons for network, lot number (22085T0005), reference (P1AU511), and serial number (007). Below this is a 'PARAMÈTRES' section with 'SIMPLE' and 'AVANCÉ' buttons. Under 'SIMPLE', there are icons for volume, labels, and signal. The 'GÉNÉRAL' section contains the following options: 'Langues' (selected, indicated by an orange arrow pointing to it), 'Tutoriel', 'IFU', 'À propos', and 'Quitter l'application'. The 'Langues' option has a sub-menu open on the right, titled 'LANGUES'. It displays the current language as 'Français' and a list of available languages: Deutsch, English, Español, Français, and Português. An 'Appliquer' button is located at the bottom right of this sub-menu.

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